

HOLOGIC®

INSTRUCTIONS FOR USE

Cleavase® X Enzyme Invader® Reagent Core Kit for Amplified DNA

Arm-Dye Code 1=FAM, 3=RED

Catalog No. 91-232

PRINCIPLES OF THE INVADER DNA ASSAY:

The Invader DNA Assays use Cleavase enzymes to recognize and cleave specific structures formed by the addition of two oligonucleotides to a nucleic acid target (Figure 1).

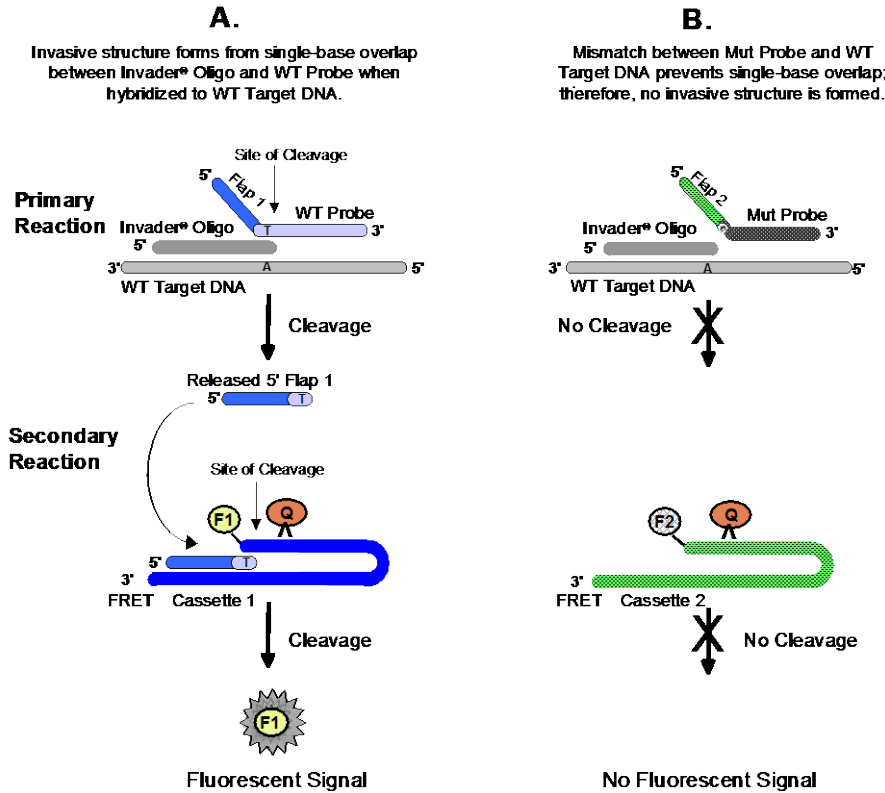


Figure 1.

Invader Reaction:

In the Invader DNA Assay, two oligonucleotides (a discriminatory Primary Probe and an Invader Oligo) hybridize in tandem to the target DNA to form an overlapping structure. The 5'-end of the Primary Probe includes a 5'-flap (arm) that does not hybridize to the target DNA (Figure 1). The 3'-nucleotide of the bound Invader Oligo overlaps the Primary Probe, but need not hybridize to the target DNA. The Cleavase enzyme recognizes this overlapping structure and cleaves off the unpaired 5'-flap (arm) of the Primary Probe, releasing it as a target-specific product. The Primary Probe is designed to have a melting temperature close to the reaction temperature. Thus, under the isothermal assay conditions, Primary Probes, which are provided in excess, cycle on the target DNA. This allows for multiple rounds of Primary Probe cleavage for each target DNA, and amplification of the number of released 5'-flaps (arms).

In the secondary reaction, each released 5'-flap (arm) can serve as an Invader Oligo on a fluorescence resonance energy transfer FRET Cassette to create another overlapping structure that is recognized and cleaved by the Cleavase enzyme (Figure 1). When the FRET Cassette is cleaved, the fluorophore (F) and quencher (Q) are separated, generating detectable fluorescence signal. Similar to the initial reaction, the released 5'-flap (arm) and the FRET Cassette cycle, resulting in amplified fluorescence signal. The initial and secondary reactions run concurrently in the same well.

The biplex format of the Invader DNA Assay enables simultaneous detection of two DNA sequences in a single well. Most often, this involves detection of two variants of a particular polymorphism. The biplex format uses two different discriminatory Primary Probes, each with a unique 5'-flap (arm), and two different FRET Cassettes, each with a spectrally distinct fluorophore. By design, the released 5'-flaps (arms) will bind only to their respective FRET Cassettes to generate a target-specific signal.

REAGENTS PROVIDED

Note: For information on any hazard and precautionary statements that may be associated with reagents, refer to the Safety Data Sheet Library at www.hologic.com/sds.

FRET Mix (2 tubes, 1.1 ml each). Dye Code Arm 1=FAM, Arm 3=RED
 Cleavase X Enzyme - MgCl₂ Mix (1 tube, 500 µl)
 No Target Blank (1 tube, 180 µl)

REAGENTS REQUIRED, BUT ORDERED SEPARATELY

Probe Mix (1.75 ml) Assay-specific probes and Invader oligonucleotide to detect the two alleles of the polymorphism of interest.

STORAGE REQUIREMENTS

Component	Storage Requirement
Probe Mix	-20°C
FRET Mix	-20°C
Cleavase X Enzyme - MgCl ₂ solution	-20°C
No Target Blank	-20°C

SUPPLIES AND EQUIPMENT REQUIRED, BUT NOT PROVIDED**Reagents:**

Mineral oil for reaction overlay to prevent evaporation. Mineral oil produced by Sigma is recommended (catalog number M-3516 or M-5904). It is imperative that a high grade of mineral oil be used, because lower grades may interfere with fluorescence readings.

Consumable Supplies:

Sterile, DNase/RNase-free disposable aerosol barrier pipet tips
 96-well or 384-well microplates
 Sterile 0.6 ml or 1.7 ml microcentrifuge tubes
 Sterile 8-tube strip or microplate (optional)
 Disposable plastic trough (optional)
 Plate sealing tape (optional)

Equipment:

Thermal cycler **OR** oven capable of holding 63° ± 1°C.
 Fluorescence plate reader (top reading; set for the appropriate excitation and emission wavelengths).
 Miscellaneous laboratory equipment (tube racks, micropipetors, multichannel pipet, microcentrifuge, vortex mixer, plate centrifuge-for manual setup only).

PRECAUTIONS

1. Use only calibrated pipets, thermal cyclers, or ovens. Do not use heat blocks, because the microtiter plates will warp, producing unreliable results.
2. Use sterile disposable aerosol barrier pipet tips for each addition and transfer to avoid cross-contamination.
3. To avoid contamination do not let the tips touch previously dispensed reagents.
4. Be sure the samples and reagents are mixed by vortexing thoroughly and consistently before use.
5. Follow the order of addition of reagents given in the **Invader DNA Assay Procedure** below.
6. Visually check the level (volume) of reagents in each plate well to verify that reagents and/or sample have been added as required.
7. Do not pool reagents from different lots or from different vials of the same lot.
8. Do not use a reagent after its expiration date.
9. Use freshly dispensed mineral oil for reaction overlay.
10. Do not transfer any reagents back to the original container once they have been dispensed.
11. Dispose of unused reagents and waste in accordance with applicable country, federal, state, and local regulations.
12. Before reading a plate, allow the lamp in the fluorescence plate reader to warm up for at least 10 minutes and verify that the plate is aligned properly in the fluorescence plate reader.
13. Take the appropriate measures necessary to prevent contamination of Invader Assay reagents, other reagents and supplies from amplified product carryover contamination.

DNA AMPLICON REQUIREMENTS

1. PCR primers should be chosen so that the amplified fragment contains at least 60 base pairs on both sides of the targeted polymorphism. Fragments of 200 to 500 base pairs in length are recommended, although longer fragments may be successfully used if the yield after amplification is adequate (see Item 5 below).
2. The PCR product must be incubated at 99°C for 10 minutes in order to inactivate the taq polymerase, which inhibits the Invader reaction. For ease of use, this step may be added at the end of the PCR program.
3. Amplify samples using optimal PCR conditions, preferably from genomic DNA. The amplified DNA may be checked on an agarose gel for integrity and yield.
4. An aliquot of the amplification reaction may be diluted and used directly in the Invader reaction. No desalting or further purification is necessary.
5. The Invader reaction will perform optimally when the PCR reaction yield is 5 to 30 fmoles/ μ l*. **PCR reactions must be diluted 1:10 to 1:50 in nuclease-free water prior to use in the Invader reaction.** After dilution, samples should contain approximately 1 to 5 fmoles/ μ l of the amplified DNA. Use of non-diluted PCR reactions is not recommended, since these inhibit the Invader reaction.

*For a 500-bp fragment, 1.6 ng/ μ l would be equivalent to 5 fmoles/ μ l.

6. The above recommendations apply to monoplex or low-plex (5-plex or less) amplification reactions.

THE INVADER DNA ASSAY PROCEDURE – 96-WELL REACTION FORMAT

The following procedure applies to assays setup manually or with automated liquid handlers. For a specific automated liquid handler system, follow the manufacturer’s recommendations.

1. Plan the microtiter plate assay layout. It is recommended to have duplicate No Target Blanks (see example below) for each assay.
2. Experimental Plate Layout: *(fill in as appropriate and required)*

	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												

An Example Layout Follows:

	1	2	3	4	5	6	7	8	9	10	11	12
A	Allele 1 Control											
B	Allele 2 Control											
C	NTB											
D	NTB											
E												
F												
G												
H												

3. Calculate and prepare the appropriate amount of Master Mix, as shown in the table below. We recommend the preparation of a master mix sufficient for testing of the complete set of samples.

Note: It is recommended to prepare a fresh mix each time.

Master Mix Component	(A) Volume per reaction	(B) Number of reactions	(C) Multiplication factor *	(A) x (B) x (C) = Volume of reagent required
Probe Mix	3 µl		1.2x	
FRET Mix	3.5 µl		1.2x	
Cleavase X Enzyme/MgCl ₂ Solution	1 µl		1.2x	
Total volume	7.5 µl		1.2x	

* This multiplication factor allows for 20% overage. For automated liquid handlers, the amount of overage may vary, depending on the dead volume and overdraws. The actual number of reactions obtained from a tube will vary, based on the liquid handler used and the number of times each tube is used.

4. Add 7.5 µl of Master Mix to each well.
5. Add 7.5 µl of the appropriate control(s) or diluted PCR product to the appropriate wells. It is recommended that the PCR product be diluted to 1 to 5 fmoles/µl in nuclease-free water (see Item 4 in DNA AMPLICON REQUIREMENTS section) which will result in a final concentration of 7.5 to 40 fmoles in the Invader reaction.

Note: Synthetic DNA Controls supplied by Hologic, Inc. (if purchased) are Invader reaction controls only. These materials are not designed to be amplification controls and therefore will not amplify.

6. Overlay each reaction with 15 µl of mineral oil.
7. Incubate the plate(s) for 20 minutes at 63°C in a thermal cycler or incubator such as the BioOven III.

Note: The incubation time at 63°C may be extended up to 60 minutes if the signal is low. Do not exceed 60 minutes.

8. Analyze the microtiter plate(s) on a fluorescence plate reader using appropriate parameters.

THE INVADER DNA ASSAY PROCEDURE – 384-WELL REACTION FORMAT

The following procedure applies to assays setup manually or with automated liquid handlers. For a specific automated liquid handler system, follow the manufacturer’s recommendations.

1. Plan the microtiter plate assay layout. It is recommended to have duplicate No Target Blanks for each assay.
2. Experimental Plate Layout: *(fill in as appropriate and required)*

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24
A																								
B																								
C																								
D																								
E																								
F																								
G																								
H																								
I																								
J																								
K																								
L																								
M																								
N																								
O																								
P																								

3. Calculate and prepare the appropriate amount of Master Mix, as shown in the table below. We recommend the preparation of a master mix sufficient for testing of the complete set of samples.

Note: It is recommended to prepare a fresh mix each time.

Master Mix Component	(A) Volume per reaction	(B) Number of reactions	(C) Multiplication factor *	(A) x (B) x (C) = Volume of reagent required
Probe Mix	1.2 µl		1.2x	
FRET Mix	1.4 µl		1.2x	
Cleavase X Enzyme/MgCl ₂ Solution	0.4 µl		1.2x	
Total volume	3 µl		1.2x	

* This multiplication factor allows for 20% overage. For automated liquid handlers, the amount of overage may vary, depending on the dead volume and overdraws. The actual number of reactions obtained from a tube will vary, based on the liquid handler used and the number of times each tube is used.

4. Add 3 µl of Master Mix to each well.
5. Add 3 µl of the appropriate control(s) or diluted PCR product to the appropriate wells. It is recommended that the PCR product be diluted to 1 to 5 fmoles/µl in nuclease-free water (see Item 4 in DNA AMPLICON REQUIREMENTS section) which will result in a final concentration of 3 to 15 fmoles in the Invader reaction.

Note: Synthetic DNA Controls supplied by Hologic, Inc. (if purchased) are Invader reaction controls only. These materials are not designed to be amplification controls and therefore will not amplify.

6. Overlay each reaction with 6.5 µl of mineral oil.
7. Incubate the plate(s) for 20 minutes at 63°C in a thermal cycler or incubator such as the BioOven III.

Note: The incubation time at 63°C may be extended up to 60 minutes if the signal is low. Do not exceed 60 minutes.

8. Analyze the microtiter plate(s) on a fluorescence plate reader using appropriate parameters listed.

DETERMINE WHICH ALLELE IS REPORTED BY FAM OR RED

Referring to the arm-allele combination listed on the probe mix vial for each assay, and the arm-dye code of the Invader Reagent Core Kit, determine which allele is reported by FAM or RED dyes:

For example:

If, the probe mix vial label lists 1=C 3=T, and the arm-dye code is 1=FAM, 3=RED
Then, the FAM dye will report the C allele, and the RED dye will report the T allele.

CALCULATION OF FOLD OVER ZERO AND ALLELIC RATIOS

1. Determine Fold Over Zero (FOZ) values as follows:

F Probe FOZ = (raw counts from F Probe for the sample) / (average raw counts from F Probe for the No Target Blank)

R Probe FOZ = (raw counts from R Probe for the sample) / (average raw counts from R Probe for the No Target Blank)

2. Calculate the Net FOZ for each sample by subtracting 1 from the FOZ values calculated in step 1 above.

Net FOZ for F Probe = F Probe FOZ -1

Net FOZ for R Probe = R Probe FOZ -1

3. Calculate the Allelic Ratio for each sample using the two Net FOZ values.

$$\text{Allelic Ratio} = \frac{\text{Net FOZ for F Probe}}{\text{Net FOZ for R Probe}}$$

NOTE: In cases where the Net FOZ value is ≤ 0 , set the value to 0.01 for the calculation of the Allelic Ratio.

RECOMMENDED INTERPRETATION CRITERIA FOR DETERMINATION OF SAMPLE GENOTYPES

Using the FOZ and the Allelic Ratio calculations, the results for each sample can be classified as:

- Homozygous for F Allele;
- Homozygous for R Allele;
- Heterozygous;
- Equivocal; or
- Invalid.

Sample determinations are valid under the following two criteria. Samples not meeting both criteria are invalid and must be repeated.

1. The Allelic Ratio, (as calculated above), falls into the ranges in the chart below and
2. The appropriate FOZ values, for a given sample, meet the criteria identified in the chart below:

Allelic Ratio	FAM FOZ (F Allele)	Redmond Red FOZ (R Allele)	Interpretation of Results
≥ 4.0	≥ 1.6	≥ 0.01	Homozygous for F Allele
≤ 0.25	≥ 0.01	≥ 1.6	Homozygous for R Allele
≥ 0.4 to ≤ 2.5	≥ 1.3	≥ 1.3	Heterozygous
> 2.5 to < 4.0	(Invalid due to Allelic Ratio)	(Invalid due to Allelic Ratio)	Equivocal (EQ1)
> 0.25 to < 0.4	(Invalid due to Allelic Ratio)	(Invalid due to Allelic Ratio)	Equivocal (EQ2)

Each laboratory may wish to establish its own ranges for Allelic Ratios for interpreting results. When equivocal results are obtained, the sample should be re-tested in duplicate. If the sample still generates an equivocal result, the concentration of the DNA sample should be verified.

Alternative methods for data analysis can be used.

TROUBLESHOOTING GUIDE

Problem	Possible Solution
No Signal or Low FOZ Signal	<p>Sample Preparation:</p> <ul style="list-style-type: none"> Insufficient amplified product was used in the assay. The amplified product sample should be quantitated using the PicoGreen® assay. Only samples with a concentration of at least 1 fmole/μl should be used in the Invader DNA Assay. The amplified product sample(s) may not have been completely denatured. Verify that the sample was denatured at the correct temperature and for an appropriate time. The DNA sample preparation method used caused amplification inhibition. A low A_{260}/A_{280} ratio (< 1.5) indicates that the sample contains excess protein. Poor amplification due to DNA degradation. Check DNA on an agarose gel to ensure that it is not degraded. Poor amplification efficiency. Verify that the amplification reaction is producing sufficient amplified product. Wrong region of the target sequence being amplified. Ensure that the PCR primers are proper for the sequence being detected. PCR primers should be chosen so that the amplified fragment contains at least 60 base pairs on both sides of the targeted polymorphism. <p>Assay:</p> <ul style="list-style-type: none"> The plate may not be oriented properly. Orient the plate so that well A-1 is in the upper left-hand corner. The wrong Invader Reagent core kit is being used with the Probe Mix. Check to see that the proper Invader Reagent core kit is being used with the Probe Mix. Mixing inconsistencies. Be sure the samples, reagents and Master Mix are mixed thoroughly and consistently. Make sure that all liquid is expelled from the tip during additions. Dispensing accuracy. For the manual format, be sure the plates have been spun down after the addition of Master Mix. For the automated format, be sure that the tip height is properly adjusted, so that solutions are delivered close to the bottom of the wells. To avoid contamination do not let the tips touch previously dispensed reagents. Verify that the correct volume of each reagent was added and that the reagents were added in the correct sequence. Amplified product and/or the Master Mix may not have been added. Verify that your pipet is functioning correctly. For the automated format, verify that you have sufficient volume in the reservoir for the automated liquid handler. The correct volume of mineral oil may not have been added. Verify that the correct volume of mineral oil was added. Lower volumes may result in sample evaporation. High variation in the No Target Blank duplicates may lead to low FOZ signal. Repeat run. The reagents were not stored properly. Verify that all reagents are stored at the proper temperature, as indicated in this package insert. Avoid subjecting reagents to repeated freeze thaw cycles. No Target Blank well(s) was/were contaminated with sample or exogenous nucleases. Take appropriate steps to avoid cross-contamination. The pipet tips may have been contaminated. Make sure that pipet tips do not touch any surfaces other than the solution being pipetted, because nucleases may be present. Do not touch pipet tips with hands. The plate may not have been incubated at 63°C for 20 – 60 minutes. Confirm that the plate was incubated properly. Poor amplification efficiency. Poor efficiency can be compensated by extending the Invader DNA Assay incubation time up to 60 minutes, in order to obtain sufficient signal. Do not exceed 60 minutes. <p>Instrument:</p> <ul style="list-style-type: none"> The plate was not oriented properly during assay setup and/or plate reading. Orient the plate so that well A-1 is in the upper left-hand corner. The gain settings were not properly adjusted. Adjust the gain of the fluorescence plate reader so that you are in the linear dynamic range of the respective reader. The lamp was not turned on long enough before reading the plate. Allow the lamp in the fluorescence plate reader to warm up for at least 10 minutes before reading the plate. The incorrect plate coordinates were used when the plate was read. Make sure the correct plate coordinates have been used, for the particular plate type being used. The probe should be set to read in the middle of the well and at an optimal distance from the plate. Consult the instrument user manual for detailed instructions or contact Hologic Technical Services for assistance. The temperature was not maintained at $63^{\circ} \pm 1^{\circ}\text{C}$. Incubations should be conducted in properly calibrated heating units. Checking these units on a regular basis using a thermocouple thermometer equipped with a probe traceable to NIST standards is recommended. Make sure that the plate is firmly seated in the thermal cycler. The plate was read at the incorrect excitation/emission wavelengths. Verify that the correct filters are installed or settings programmed properly for the respective fluorescence plate reader being used. If possible, run a diagnostic test on the fluorescence plate reader to ensure that the instrument and light source are working properly. Verify that two scans were performed at the recommended wavelengths.

TROUBLESHOOTING GUIDE - CONTINUED

Problem	Possible Solution
High Signal in the No Target Blank	<p>Assay:</p> <ul style="list-style-type: none"> • Use DNase/RNase free aerosol barrier tips and sterile tubes for making the Master Mix. • Wear gloves when setting up the assay. • Make sure that pipet tips do not touch any other surfaces except the solution being pipetted, since nucleases may be present. • Do not touch pipet tips with hands. • No Target Blank may be contaminated with amplified product. Use an unopened vial of No Target Blank and repeat the assay.
Fluorescent Signal Is Off-scale	<p>Assay:</p> <ul style="list-style-type: none"> • Use DNase/RNase free aerosol barrier tips and sterile tubes for making the Master Mix. • Confirm that the incubations were done for the correct amount of time and at the correct temperature. • Amplified product was not diluted sufficiently. Confirm that amplification reaction was diluted 1:10 to 1:50 to yield a concentration of approximately 1 to 5 fmoles/μl of the amplified DNA. • The Invader reaction incubation time was too long. Do not incubate the reactions for longer than 60 minutes. • The wrong Invader Reagent Core Kit is being used. Ensure that the Invader Reagent Core Kit for Amplified DNA is being used. <p>Instrument:</p> <ul style="list-style-type: none"> • Adjust the gain of the fluorescence plate reader so that you are in the linear dynamic range of the respective reader. • Allow the lamp in the fluorescence plate reader to warm up for at least 10 minutes before reading the results.
The No Target Blank Values for the F dye and R dye Are Not the Same.	<p>Instrument:</p> <ul style="list-style-type: none"> • The No Target Blank values for the F and R dye scans will vary based on the gain setting and from assay to assay. It is not uncommon for the No Target Blank values to be different for the F and R dye scans.
The Net FOZ Value from Samples Is Less Than 0.	<p>Assay:</p> <ul style="list-style-type: none"> • Samples can produce a Net FOZ value less than 0 for the non-reactive probe. This may be due to differences in the matrix of the sample compared to the matrix of the No Target Blank. NOTE: In cases where the Net FOZ value is ≤ 0, set the value to 0.01 for the calculation of the Allelic Ratio.
The Allelic Ratio of the Sample Is Close to the Limit of the Recommended Value or Equivocal.	<p>Assay:</p> <ul style="list-style-type: none"> • The Allelic Ratios provided should only be used as a guideline. • High number of equivocal calls due to non-specific target amplification. Verify primer design to ensure specific target amplification. • High number of equivocal calls due to non-specific signal. Amplified product was not diluted sufficiently. Confirm that amplification reaction was diluted 1:10 to 1:50 to yield a concentration of approximately 1 to 5 fmoles/μl of the amplified DNA. • Verify primer design to ensure specific target amplification. • Check the concentration of DNA using the PicoGreen assay since low FOZ can lead to ratios near the limits or to Equivocal results.

REFERENCES:

1. Nagano, M. *et al.*, Two novel missense mutations in the CETP gene in Japanese hyperalphalipoproteinemic subjects: high-throughput assay by Invader® assay. *Journal of Lipid Research* 43, 1011-1018 (2002)
2. Neville, M. *et al.*, Characterization of cytochrome P450 alleles using the Invader System. *Biotechniques* 32, S34-S43 (2002)
3. Hessner, M. J. *et al.*, Multisite study for genotyping of the factor II (prothrombin) G20210A mutation by the Invader assay. *Clin Chem.* 47, 2048-2050 (2001).
4. Hall, J. G. *et al.*, Sensitive detection of DNA polymorphisms by the serial invasive signal amplification reaction. *Proc Natl Acad Sci U S A.* 97, 8272-8277 (2000).
5. Hessner, M. J. *et al.*, Genotyping of factor V G1691A (Leiden) without the use of PCR by invasive cleavage of oligonucleotide probes. *Clin Chem.* 46, 1051-1056 (2000).
6. Ledford, M. *et al.*, A multi-site study for detection of the factor V (Leiden) mutation from genomic DNA using a homogeneous Invader® microtiter plate FRET assay. *J Molec Diagnostics.* 2, 97-104 (2000).
7. Lyamichev, V. I. *et al.*, Experimental and theoretical analysis of the invasive signal amplification reaction. *Biochemistry.* 39, 9523-9532 (2000).
8. Reynaldo, L. P. *et al.*, The kinetics of oligonucleotide replacements. *J Mol Biol.* 297, 511-520 (2000).
9. Kaiser, M. W. *et al.*, A comparison of eubacterial and archaeal structure-specific 5'- exonucleases. *J Biol Chem.* 274, 21387-21394 (1999).
10. Kwiatkowski, R. W. *et al.*, Clinical, genetic, and pharmacogenetic applications of the Invader™ assay. *Mol Diagn.* 4, 353-364 (1999).
11. Lyamichev, V. *et al.*, Polymorphism identification and quantitative detection of genomic DNA by invasive cleavage of oligonucleotide probes. *Nat Biotechnol.* 17, 292-296 (1999).
12. Patnaik M. *et al.* Detection of genomic polymorphisms associated with venous thrombosis using the Invader biplex assay. *J Mol Diagn.* 6, 137-144 (2004).

Manufactured and Distributed by Hologic, Inc., San Diego, CA USA

For further technical information, contact:

Hologic, Inc.
10210 Genetic Center Drive
San Diego, CA 92121 USA

Customer Support: +1 844 Hologic (+1 844 465 6442)
customersupport@hologic.com

Technical Support: +1 888 484 4747
molecularsupport@hologic.com

For more contact information visit www.hologic.com.

NOTICE TO RECIPIENT ABOUT LIMITED LICENSE

The receipt of this product from Hologic, Inc. or its authorized distributor includes a limited, non-exclusive license under patent rights held by Hologic, Inc. Acquisition of this product constitutes acceptance by the recipient of this limited license. Recipients unwilling to accept the limited license must return the product for a full refund. Such license is solely for the purposes of using this product to detect a specific analyte. For avoidance of doubt, the foregoing license does not include rights to use this product for agriculture or veterinary medicine applications. The foregoing license does not include a license to use the product for new product research or development, product manufacture, or any reverse-engineering purposes. The purchaser of this product is not authorized to transfer this product to any third party for any purpose without the express written consent of Hologic, Inc. Except as expressly provided in this paragraph, no other license is granted expressly, impliedly, or by estoppel. For information concerning the availability of additional licenses to practice the patented methodologies, contact: Legal Department, Hologic, Inc., 250 Campus Drive, Marlborough, MA, 01752, (508) 263-2900.

This product may be covered by one or more U.S. patents identified at www.hologic.com/patents.

LIMITED PRODUCT WARRANTY

WARRANTIES. Equipment, Supplies, and Software are warranted to the original Customer to perform substantially in accordance with published Product Specifications for one (1) year starting from the date of Installation (if applicable) or from the date of Delivery, whichever occurs first. After-sale options and accessories are warranted for six (6) months, and x-ray tubes are warranted on a straight-line prorated basis as stated in the applicable Product Specification ("Warranty Period"). Replacement parts are warranted for the remainder of the Warranty Period or ninety (90) days from Delivery, whichever is longer. Consumable Supplies are warranted to conform to published specifications for a period ending on the expiration date shown on their respective packages. Services are warranted to be supplied in a workman-like manner. Hologic does not warrant that use of Products will be uninterrupted or error-free, or that Products will operate with non-Hologic authorized third-party products. HOLOGIC'S ENTIRE WARRANTY RESPONSIBILITY IS EXPRESSLY LIMITED TO REPAIR OR REPLACEMENT (AT HOLOGIC'S OPTION AND IN THE FORM ORIGINALLY SHIPPED) OF PRODUCT OR CORRECTION OF SERVICE SUBJECT TO ANY CLAIM, OR, AT HOLOGIC'S ELECTION, REPAYMENT OF, OR CREDITING CUSTOMER WITH, AN AMOUNT EQUAL TO THE HOLOGIC PRICE, FEE OR CHARGE THEREFORE. THE FOREGOING WARRANTIES ARE IN LIEU OF AND EXCLUDE ALL OTHER WARRANTIES NOT EXPRESSLY SET FORTH HEREIN, WHETHER EXPRESS OR IMPLIED BY OPERATION OF LAW OR OTHERWISE, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. SUCH LIMITED WARRANTY IS GIVEN SOLELY TO THE ORIGINAL CUSTOMER AND IS NOT GIVEN TO, NOR MAY IT BE RELIED UPON BY, ANY THIRD PARTY INCLUDING, WITHOUT LIMITATION, CUSTOMERS OF CUSTOMER. THIS WARRANTY IS VOID UPON TRANSFER OF PRODUCT BY CUSTOMER TO ANY ENTITY WHO HAS LESS THAN FIFTY (50) PERCENT OWNERSHIP IN THE PRODUCT. SOME STATES DO NOT ALLOW THE EXCLUSION OF IMPLIED WARRANTIES SO THE ABOVE EXCLUSIONS MAY NOT APPLY TO YOU. YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY, FROM STATE TO STATE. These warranties do not apply to any item that is: (a) repaired, moved or altered other than by Hologic authorized service personnel; (b) subjected to physical (including thermal or electrical) abuse, stress, or misuse; (c) stored, maintained or operated in any manner inconsistent with applicable Hologic specifications or instructions; or (d) designated as supplied subject to a non-Hologic warranty or on a pre-release or "as-is" basis.

WARRANTY CLAIMS AND REMEDIES. In the event of any warranty claim, Hologic will replace with new or repaired items any Equipment part, component, or consumable supply that is in breach of warranty, and will use reasonable efforts to promptly fix or provide a workaround for any Software defect or bug which prevents operation in substantial conformity with functional specifications. Alternatively, Hologic may elect to repay or credit to Customer an amount equal to the purchase price of the defective Equipment, component, Software, consumable supply, or Service. Items replaced shall become Hologic property. All claims shall be initiated by contacting Hologic within the applicable warranty period and thirty (30) days after discovery of the breach or non-conformity. Hologic must be given reasonable access and an opportunity to inspect all associated materials. If Hologic and Customer are unable to settle any claim and Customer has not notified Hologic within one (1) year after the claim arises, Customer shall be barred from instituting any legal action thereafter. These remedies shall comprise Hologic's entire liability and Customer's exclusive remedy for breach of warranty and are in lieu of any other remedies at law or equity.

LIMIT OF LIABILITY. HOLOGIC SHALL NOT BE LIABLE FOR ANY SPECIAL, INCIDENTAL, PUNITIVE, EXEMPLARY OR CONSEQUENTIAL LOSSES, DAMAGES, OR EXPENSES (INCLUDING BUT NOT LIMITED TO LOSS OF PROFITS, DATA, OR USE), DIRECTLY OR INDIRECTLY ARISING FROM THE SALE, HANDLING, SERVICE OR USE OF PRODUCT ORDERED OR FURNISHED, OR FROM ANY CAUSE RELATING THERETO UNLESS EXPRESSLY AGREED TO BY THE PARTIES IN WRITING. EXCEPT FOR PERSONAL INJURY OR DEATH TO THE EXTENT RESULTING FROM HOLOGIC'S NEGLIGENT OR INTENTIONALLY WRONGFUL ACTS OR OMISSIONS, IN NO EVENT SHALL HOLOGIC BE LIABLE UNDER ANY LEGAL THEORY OR FOR ANY CAUSE WHATSOEVER, WHETHER BASED UPON WARRANTY, CONTRACT, TORT, NEGLIGENCE, OR OTHER THEORY, EVEN IF ADVISED OF THE POSSIBILITY THEREOF, FOR ANY AMOUNT IN EXCESS OF THE PRICE, FEE OR CHARGE THEREFORE RECEIVED BY HOLOGIC

Hologic, Cleavase, and Invader are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries.

All other trademarks that may appear in this package insert are the property of their respective owners.

Some components of nucleic acid analysis, such as specific methods and compositions for manipulating or visualizing nucleic acids for analysis, may be covered by one or more patents owned by other parties. Similarly, nucleic acids containing specific nucleotide sequences may be patented. Making, using, offering for sale, or selling such components or nucleic acids may require one or more licenses. Nothing in this document should be construed as an authorization or implicit license to make, use or sell any so covered component or nucleic acid under any such patents.

©2011-2015 Hologic, Inc. All rights reserved.
Part Number 15-1592, revision 101