

Hologic® RespDirect™ Collection Kit

For *in vitro* diagnostic use.

Rx only

Intended Use

The Hologic® RespDirect™ Collection Kit is intended to be used for the collection of nasopharyngeal (NP) swab specimens (collected by a healthcare professional) for testing with the Panther Fusion SARS-CoV-2/Flu A/B/RSV assay.

Materials Provided

50 Hologic RespDirect Collection Kits (Cat. No. PRD-07788)

Each kit contains:

Component	Quantity	Description
RespDirect Swab	1	Individually wrapped, sterile swab
Enhanced Direct Load Tube (eDLT)	1	Tube containing Enhanced Specimen Transport Medium (eSTM), 2.9 mL

Warnings and Precautions

- A. Use the provided swab only. Failure to use the provided swab may invalidate the test results.
- B. Do not use if the swab, tube, or packaging is visibly damaged (e.g., if the swab tip or shaft is bent or broken).
- C. Do not bend or shape the swab before collection. Do not use excessive force, pressure, or bending when collecting the specimen as this may result in accidental breakage of the swab.
- D. Do not apply the eSTM directly to skin or mucous membranes or take internally.
- E. Specimens may be infectious. Use Universal Precautions when handling specimens. Only personnel adequately trained in handling infectious materials should be permitted to handle specimens.
- F. Take care to avoid cross-contamination during specimen handling steps. Specimens can contain extremely high levels of pathogens. Ensure that specimens do not contact one another, and discard used materials without passing over the specimen tubes. If gloves come into contact with a specimen, change gloves to avoid cross-contamination.
- G. If the contents of the tube are spilled at any time during the collection procedure, collect a new specimen using a new Hologic RespDirect Collection Kit. Failure to use a new kit may invalidate the test results.
- H. Maintain proper storage conditions during specimen shipping to ensure the integrity of the specimen. Specimen stability under shipping conditions other than those recommended has not been evaluated.
- I. Do not use the kit after the expiration date to collect specimens.

- J. Use appropriate personal protective equipment (PPE) when collecting and handling specimens from individuals suspected of being infected with SARS-CoV-2 as outlined in CDC Interim Laboratory Biosafety Guidelines for Handling and Processing Specimens Associated with Coronavirus Disease 2019 (COVID-19).¹
- K. Avoid contamination by adhering to good laboratory practices and to the procedures specified in this package insert.
- L. Some reagents in the kit are labeled with hazard information.

Note: Hazard communication information for labeling of globally marketed products reflects the US Safety Data Sheets (SDS) classifications. For hazard communication information specific to your region, refer to the region specific SDS on the Safety Data Sheet Library at www.hologicds.com. For more information on the symbols, refer to the symbol legend on www.hologic.com/package-inserts.

US Hazard Information	
	<p>Enhanced Specimen Transport Media HEPES 1-5% Subtilisins (proteolytic enzymes) <1% 1,2-Benzisothiazolin-3-one <1%</p> <p>DANGER May cause allergy or asthma symptoms or breathing difficulties if inhaled Avoid breathing dust/fume/gas/mist/vapors/spray In case of inadequate ventilation wear respiratory protection Wear respiratory protection In case of inadequate ventilation wear respiratory protection IF INHALED: Remove person to fresh air and keep comfortable for breathing IF INHALED: Remove person to fresh air and keep comfortable for breathing IF INHALED: Remove person to fresh air and keep comfortable for breathing IF INHALED: Remove person to fresh air and keep comfortable for breathing</p>

Limitations

- A. The Hologic RespDirect Collection Kit is not cleared/indicated for home collection.

Kit Storage Requirements

Prior to use store collection kit at room temperature (15°C to 30°C).

RespDirect Specimen Performance

Assay performance characteristics for specimens collected with the RespDirect collection kit are provided in the Panther Fusion SARS-CoV-2/Flu A/B/RSV Assay package insert.


Nasopharyngeal Swab Specimens Collection and Handling

Instructions for nasopharyngeal swab specimen collection (healthcare professional collection only):

1. Wash hands with soap and water. If soap and water are not available, use hand sanitizer.
2. Partially peel open the swab package. Remove the swab. Do not touch the soft tip or lay the swab down. If the soft tip is touched, laid down, or dropped, discard and use a new Hologic RespDirect Collection Kit.
3. Hold the swab, placing your thumb and forefinger below the score line.


¹ [https://www.who.int/publications/i/item/laboratory-biosafety-guidance-related-to-coronavirus-disease-\(covid-19\)](https://www.who.int/publications/i/item/laboratory-biosafety-guidance-related-to-coronavirus-disease-(covid-19))

Note: Actual product may appear slightly different from the one shown in this section.

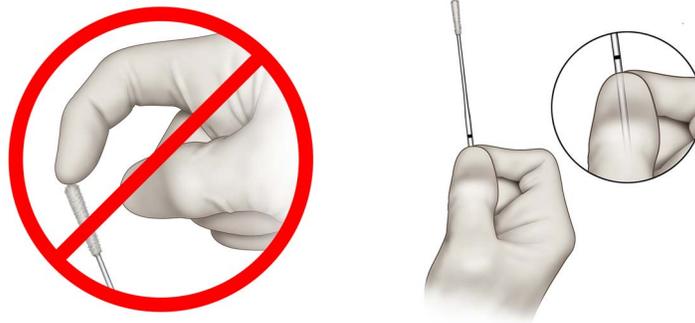

4. Tilt patient's head back 70 degrees. Gently and slowly insert the swab through the nostril parallel to the palate (not upwards) until resistance is encountered or the distance is equivalent to that from the ear to the nostril of the patient, indicating contact with the nasopharynx. Swab should reach the depth equal to the distance from nostrils to outer opening of the ear. Gently rub and roll the swab. Leave the swab in place for several seconds to absorb secretions. Slowly remove the swab while rotating it.²

Note: If the swab is saturated with fluid from the first collection, it is not necessary to collect specimens from both sides. If a deviated septum or blockage creates difficulty in obtaining the specimen from one nostril, use the same swab to obtain the specimen from the other nostril.

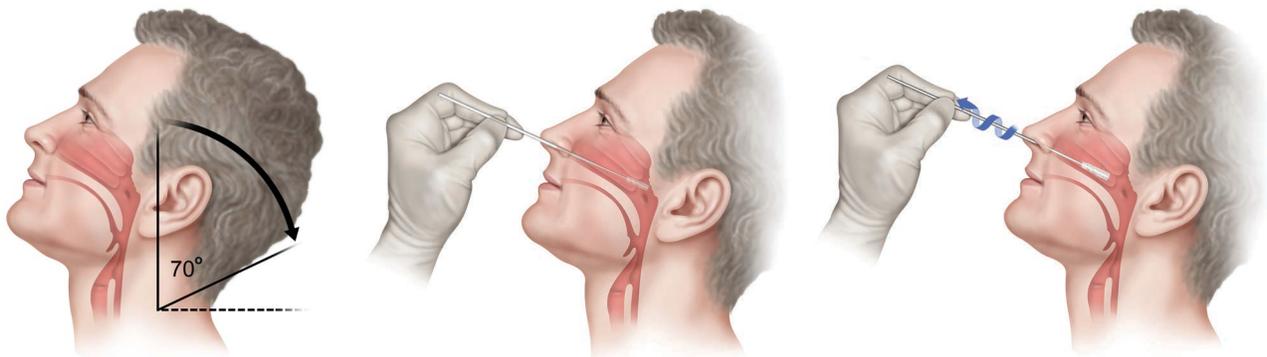


5. While holding the swab in the same hand, unscrew the cap from the tube. Do not spill the contents of the tube. If the contents of the tube are spilled, collect a new specimen using a new Hologic RespDirect Collection Kit.

Warning: If your skin comes into contact with liquid from the tube, wash with soap and water for 30 seconds.


6. Immediately place the swab into the tube so that the score line is at the top of the tube.

² <https://www.cdc.gov/coronavirus/2019-ncov/lab/guidelines-clinical-specimens.html>

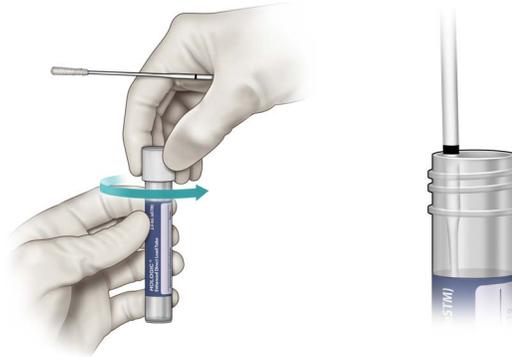

7. Carefully break the swab shaft at the score line against the side of the tube.
8. Discard the top portion of the swab shaft.
9. Tightly screw the cap onto the tube.

Note: *The tube label includes a field to record the specimen source and additional patient information.*

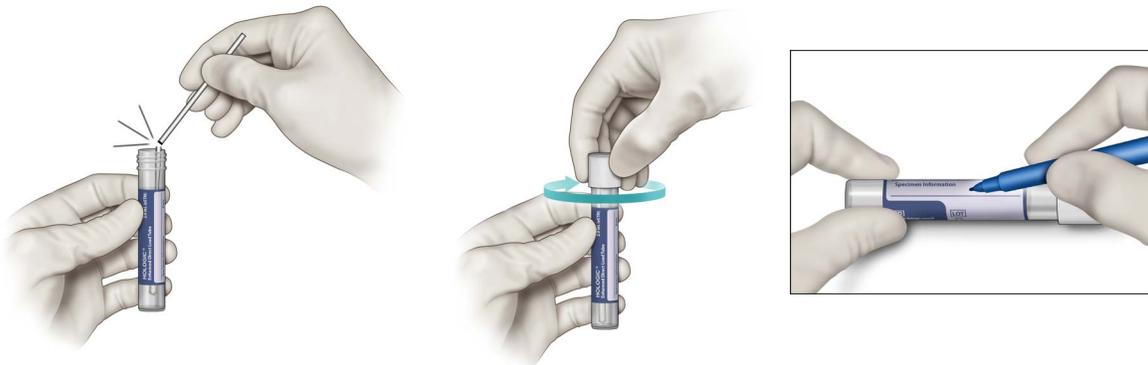

Specimen Transport and Storage

Swab specimens are intended to be transported to the laboratory in the provided swab specimen transport medium and tube. Refer to the Panther Fusion SARS-CoV-2/Flu A/B/RSV assay package insert for proper specimen storage conditions.

Note: *Specimens must be shipped in accordance with applicable national and international regulations.*

Contact Information



10210 Genetic Center Drive
San Diego, CA 92121 USA

Customer Support: +1 800 442 9892
customersupport@hologic.com

Technical Support: +1 888 484 4747
molecularsupport@hologic.com

For more contact information visit www.hologic.com.

Hologic and associated logos are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries.

All other trademarks that may appear in this package insert are the property of their respective owners.

This product may be covered by one or more U.S. patents identified at www.hologic.com/patents.

©2023 Hologic, Inc. All rights reserved.
AW-27716-001 Rev. 001
2023-05