

Aptima® HBV Quant Assay

For *in vitro* diagnostic use Rx Only

General Information	2
Intended Use	2
Summary and Explanation of the Test	2
Principles of the Procedure	3
Warnings and Precautions	3
Reagent Storage and Handling Requirements	6
Specimen Collection and Storage	7
Samples Onboard the Panther System	9
Specimen Transport	10
Panther System	11
Reagents and Materials Provided	
Materials Required but Available Separately	
Optional Materials	
Panther System Test Procedure	
Procedural Notes	
Quality Control	18
Assay Calibration	
Negative and Positive Controls	
Internal Calibrator/Internal Control	18
Interpretation of Results	19
Limitations	
Nonclinical Performance	20
Nonclinical Performance	
Limit of Detection Using the 3rd WHO International Standard	20
Limit of Detection Using the 3rd WHO International Standard	
Limit of Detection Using the 3rd WHO International Standard	
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes	
Limit of Detection Using the 3rd WHO International Standard	
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard	
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes	
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard	20 21 22 22 22 24 26
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision	20 21 22 22 22 24 26 26
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances	20 21 22 22 24 26 26 28
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity	20 20 21 22 22 24 26 26 28
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100)	20 21 22 22 24 26 26 28 29 31
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent	20 21 22 22 24 26 26 28 31 31
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent Precision of Diluted Samples	20 21 22 22 24 26 26 28 29 31 31 31 32
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent Precision of Diluted Samples Carryover	20 21 22 22 24 26 26 28 31 31 32 33 33
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent Precision of Diluted Samples	20 21 22 22 24 26 26 28 31 31 32 33 33
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent Precision of Diluted Samples Carryover Clinical Performance Clinical Utility	20 20 21 22 24 26 26 28 31 31 31 32 33 33 34 35
Limit of Detection Using the 3rd WHO International Standard Limit of Detection Across HBV Genotypes Linear Range Linearity Across HBV Genotypes Lower Limit of Quantitation Using the 3rd WHO International Standard Lower Limit of Quantitation Across HBV Genotypes Traceability to the 3rd WHO International Standard Precision Potentially Interfering Substances Analytical Specificity Matrix Equivalency Sample Dilution Using Aptima Specimen Diluent (1:3) Sample Dilution Using Aptima Specimen Diluent (1:100) Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent Precision of Diluted Samples Carryover Clinical Performance	20 21 22 22 24 26 26 28 31 31 32 33 33 34 34

General Information Aptima®

General Information

Intended Use

The Aptima HBV Quant assay is an *in vitro* nucleic acid amplification test for the quantitation of hepatitis B virus (HBV) DNA in human plasma and serum on the fully automated Panther® system.

Plasma may be prepared in ethylenediaminetetraacetic acid (EDTA), anticoagulant citrate dextrose (ACD) solution, and plasma preparation tubes (PPTs). Serum may be prepared in serum tubes and serum separator tubes (SSTs). Specimens are tested using the fully automated Panther system for sample processing, amplification, and quantitation. Specimens containing HBV genotypes A, B, C, D, E, F, G, and H are validated for quantitation in the assay.

The Aptima HBV Quant assay is intended for use as an aid in the management of patients with chronic HBV infections undergoing HBV antiviral drug therapy. The assay can be used to measure HBV DNA levels at baseline and during treatment to aid in assessing viral response to treatment. The results from the Aptima HBV Quant assay must be interpreted within the context of all relevant clinical and laboratory findings. Assay performance for determining the clinical stage of HBV infection has not been established. Clinical performance characteristics have been established for individuals treated with tenofovir disoproxil fumarate or entecavir.

The Aptima HBV Quant assay is not approved for use as a screening test for the presence of HBV DNA in blood or blood products or as a diagnostic test to confirm the presence of HBV infection.

Summary and Explanation of the Test

Hepatitis B virus, one of several viruses known to cause hepatitis, has been attributed to lifelong HBV infection, cirrhosis of the liver, liver cancer, liver failure, and potentially, death. The World Health Organization (WHO) lists HBV as one of the world's most common infectious diseases. The prevalence of HBV infection and method of transmission varies greatly around the world. About one third of the world's population has serological evidence of past or present HBV infection with chronic HBV infection occurring in an estimated 240 million people worldwide. HBV infection results in increased risk of hepatic decompensation, cirrhosis, and hepatocellular carcinoma (HCC) with a mortality of 0.5 to 1.2 million deaths and 5-10% of cases of liver transplantation worldwide annually. Without appropriate treatment, intervention, and monitoring after diagnosis, the 5 year cumulative incidence of cirrhosis ranges from 8-20%. Once cirrhosis has developed, the annual risk of hepatocellular carcinoma is 2-5%.

HBV contains a circular, partially double-stranded DNA genome of approximately 3200 base pairs, which encode four partially overlapping open reading frames (ORF) expressing the polymerase, surface, precore/core, and X proteins. The polymerase ORF overlaps the other three ORFs and encodes a key viral replication protein, polymerase. The surface ORF expresses three proteins, which are essential for viral morphogenesis, viral entry into hepatocytes, and provoking the host's immune response. There are 8 HBV genotypes (A-H), and these are typically found in distinct geographical locations.

Due to the dynamic nature of chronic hepatitis B infection, continuous monitoring of HBV DNA and alanine aminotransferase (ALT) levels is important. For the majority of HBV infected individuals who are undergoing antiviral treatment, the goal is HBV DNA suppression. Quantitative nucleic acid tests with a broad linear range are effective tools to monitor HBV DNA viral load during the course of treatment.

Aptima® General Information

Principles of the Procedure

The Aptima HBV Quant assay is an in vitro nucleic acid amplification test that uses real time transcription-mediated amplification (TMA) technology on the Panther system to quantitate HBV DNA, genotypes A, B, C, D, E, F, G, and H. The Aptima HBV Quant assay targets two highly conserved regions in the polymerase and surface genes (for increased tolerance to potential mutations). The assay is standardized to the 3rd WHO International Standard for Hepatitis B Virus (NIBSC code: 10/264).

The Aptima HBV Quant assay involves three main steps, which all take place in a single tube on the Panther system: target capture, target amplification by TMA, and detection of the amplification products (amplicon) by the fluorescent labeled probes (torches).

During target capture, viral DNA is isolated from specimens. The specimen is treated with a detergent to solubilize the viral envelope, denature proteins, and release viral genomic DNA. Capture oligonucleotides hybridize to highly conserved regions of HBV DNA, if present, in the test specimen. The hybridized target is then captured onto magnetic microparticles that are separated from the specimen in a magnetic field. Wash steps remove extraneous components from the reaction tube.

Target amplification occurs via TMA, which is a transcription-mediated nucleic acid amplification method that utilizes two enzymes, Moloney murine leukemia virus (MMLV) reverse transcriptase and T7 RNA polymerase. The reverse transcriptase is used to generate a DNA copy (containing a promoter sequence for T7 RNA polymerase) of the target sequence. T7 RNA polymerase produces multiple copies of RNA amplicon from the DNA copy template. The Aptima HBV Quant assay utilizes the TMA method to amplify two regions of the HBV genome (polymerase gene and surface gene). Amplification of those regions is achieved using specific primers designed to amplify HBV genotypes A, B, C, D, E, F, G, and H. The dual target region approach with primer design targeting the highly conserved regions ensures accurate quantitation of the HBV DNA.

Detection is achieved using single-stranded nucleic acid torches that are present during the amplification of the target, which hybridize specifically to the amplicon in real time. Each torch has a fluorophore and a quencher. When the torch is not hybridized to the amplicon, the guencher is in close proximity of the fluorophore and suppresses the fluorescence. When the torch binds to the amplicon, the quencher is moved farther away from the fluorophore and will emit a signal at a specific wavelength when excited by a light source. As more torches hybridize to amplicon, a higher fluorescent signal is generated. The time taken for the fluorescent signal to reach a specified threshold is proportional to the starting HBV concentration. Each reaction has an internal calibrator/internal control (IC) that controls for variations in specimen processing, amplification, and detection. The concentration of a sample is determined by the Panther system software using the HBV and IC signals for each reaction and comparing them to calibration information.

Warnings and Precautions

- A. To reduce the risk of invalid results, carefully read the entire package insert and the Panther System Operator's Manual prior to performing this assay.
- B. Target Enhancer Reagent (TER) is corrosive. See Safety Data Sheet information at the end of this section.

Laboratory Related



C. CAUTION: The controls for this assay contain human plasma. The plasma is non-reactive for hepatitis B surface antigen (HBsAg), antibodies to HCV, antibodies to HIV-1 and

General Information Aptima®

HIV-2, and HIV antigen when tested with US Food and Drug Administration licensed procedures. In addition, the plasma is nonreactive for HBV DNA, HCV RNA, and HIV-1 RNA when tested with licensed nucleic acid tests. All human blood sourced materials should be considered potentially infectious and should be handled with Universal Precautions.^{8,9,10}

- D. Only personnel adequately trained in the use of the Aptima HBV Quant assay and in handling potentially infectious materials should perform this procedure. If a spill occurs, immediately disinfect following appropriate site procedures.
- E. Use only supplied or specified disposable laboratory ware.
- F. Use routine laboratory precautions. Do not pipet by mouth. Do not eat, drink or smoke in designated work areas. Wear disposable, powderless gloves, protective eye wear, and laboratory coats when handling specimens and kit reagents. Wash hands thoroughly after handling specimens and kit reagents.
- G. Work surfaces, pipettes, and other equipment must be regularly decontaminated with 2.5% to 3.5% (0.35 M to 0.5 M) sodium hypochlorite solution.
- H. Dispose of all materials that have come in contact with specimens and reagents according to local, state, and federal regulations.^{8,9,10,11} Thoroughly clean and disinfect all work surfaces.
- I. The controls contain sodium azide as a preservative. Do not use metal tubing for reagent transfer. If solutions containing sodium azide compounds are disposed of in a plumbing system, they should be diluted and flushed with generous amounts of running water. These precautions are recommended to avoid accumulation of deposits in metal piping in which explosive conditions could develop.
- J. Good standard practices for molecular laboratories include environmental monitoring. To monitor a laboratory's environment, the following procedure is suggested:
 - 1. Obtain a cotton-tipped swab and pair with the Aptima Specimen Aliquot Tube (SAT).
 - 2. Label each SAT appropriately.
 - 3. Fill each SAT with 1 mL of Aptima Specimen Diluent.
 - 4. To collect the surface samples, lightly moisten a swab with nuclease-free deionized water.
 - 5. Swab the surface of interest using a top to bottom vertical motion. Rotate the swab approximately one-half turn while swabbing the location.
 - 6. Immediately place the swab sample into the tube and gently swirl the swab in the diluent to extract potential swabbed materials. Press the swab on the side of the transport tube to extract as much liquid as possible. Discard the swab and cap the tube.
 - 7. Repeat steps for remaining swab samples.
 - 8. Test swab with molecular assay.

Specimen Related

- K. Specimens may be infectious. Use Universal Precautions^{8,9,10} when performing this assay. Proper handling and disposal methods should be established in accordance with local regulations.¹¹ Only personnel adequately trained in the use of the Aptima HBV Quant assay and trained in handling infectious materials should perform this procedure.
- L. Maintain proper storage conditions during specimen shipping to ensure the integrity of the specimen. Specimen stability under shipping conditions other than those recommended has not been evaluated.

M. Avoid cross-contamination during the specimen handling steps. Be especially careful to avoid contamination by the spread of aerosols when loosening or uncapping specimens. Specimens can contain extremely high levels of organisms. Ensure that specimen containers do not contact one another, and discard used materials without passing over open containers. Change gloves if they come in contact with specimen.

Assay Related

- N. Do not use the reagent kit, the calibrator, or the controls after the expiration date.
- O. Do not interchange, mix, or combine assay reagents from kits with different master lot numbers. Assay fluids can be from different lot numbers. Controls and the calibrator can be from different lot numbers.
- P. Avoid microbial and nuclease contamination of reagents.
- Q. Cap and store all assay reagents at specified temperatures. The performance of the assay may be affected by use of improperly stored assay reagents. See *Reagent Storage* and *Handling Requirements* and *Panther System Test Procedure* for more information.
- R. Do not combine any assay reagents or fluids without specific instruction. Do not top off reagents or fluids. The Panther system verifies reagent levels.
- S. Avoid contact of TER with skin, eyes, and mucous membranes. Wash with water if contact with this reagent occurs. If spills of this reagent occur, dilute with water and follow appropriate site procedures.
- T. Some reagents of this kit are labeled with risk and safety symbols.

Note: For hazard communication information, refer to the Safety Data Sheet Library at www.hologicsds.com.



Target Enhancer Reagent

LITHIUM HYDROXIDE, MONOHYDRATE 5-10%

DANGER



H302 – Harmful if swallowed.

H314 – Causes severe skin burns and eye damage.

P264 - Wash face, hands and any exposed skin thoroughly after handling

P270 - Do not eat, drink or smoke when using this product

P260 - Do not breathe dust/fume/gas/mist/vapors/spray

P280 - Wear protective gloves/protective clothing/eye protection/face protection

P301 + P330 + P331 - IF SWALLOWED: rinse mouth. Do NOT induce vomiting

P301 + P312 - IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell

P303 + P361 + P353 – IF ON SKIN (or hair): Remove/Take off immediately all contaminated clothing. Rinse skin with water/shower

P304 + P340 - IF INHALED: Remove victim to fresh air and keep at rest in a position comfortable for breathing

P305 + P351 + P338 - IF IN EYES: Rinse cautiously with water for several minutes.

Remove contact lenses, if present and easy to do. Continue rinsing

P310 - Immediately call a POISON CENTER or doctor/physician

P330 - Rinse mouth

P363 - Wash contaminated clothing before reuse

P405 - Store locked up

Dispose of contents/container to an approved waste disposal plant.

General Information Aptima®



Controls Kit

SODIUM AZIDE 0.2%

Defibrinated Human Plasma 95-100%



WARNING

EUH032 - Contact with acids liberates very toxic gas

H302 - Harmful if swallowed

H312 - Harmful in contact with skin.

H402 - Harmful to aquatic lifeH412 - Harmful to aquatic life with long lasting effects.

P280 - Wear protective gloves/protective clothing/eye/face protection

P301 + P312 - IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell

P302 + P352 - IF ON SKIN: Wash with plenty of soap and water

P363 - Wash contaminated clothing before reuse

P501 - Dispose of contents/container in accordance with local/regional/national/international regulation.

Reagent Storage and Handling Requirements

A. The following table shows the storage conditions and stability for reagents, controls, and calibrator.

Promot	Unopened	Open Kit (Reconstituted)		
Reagent	Storage	Storage	Stability	
qHBV Amplification Reagent	2°C to 8°C	N/A	N/A	
qHBV Amplification Reconstitution Solution	2°C to 8°C	2°C to 8°C	30 days ^a	
qHBV Enzyme Reagent	2°C to 8°C	N/A	N/A	
qHBV Enzyme Reconstitution Solution	2°C to 8°C	2°C to 8°C	30 days ^a	
qHBV Promoter Reagent	2°C to 8°C	N/A	N/A	
qHBV Promoter Reconstitution Solution	2°C to 8°C	2°C to 8°C	30 days ^a	
qHBV Target Capture Reagent	2°C to 8°C	2°C to 8°C	30 days ^a	
qHBV PCAL (Positive Calibrator)	-15°C to -35°C	15°C to 30°C	Single use vial Use within 24 hours	
qHBV NC CONTROL – (Negative Control)	-15°C to -35°C	15°C to 30°C	Single use vial Use within 24 hours	
qHBV LPC CONTROL + (Low Positive Control)	-15°C to -35°C	15°C to 30°C	Single use vial Use within 24 hours	
qHBV HPC CONTROL + (High Positive Control)	-15°C to -35°C	15°C to 30°C	Single use vial Use within 24 hours	
qHBV Target Enhancer Reagent	15°C to 30°C	15°C to 30°C	30 days ^a	

^a When reagents are removed from the Panther system, they should be immediately returned to their appropriate storage temperatures.

- B. Discard any unused, reconstituted reagents, target capture reagent (TCR), and target enhancer reagent (TER) after 30 days or after the Master Lot expiration date, whichever comes first.
- C. Reagents stored onboard the Panther system have 72 hours of onboard stability. Reagents can be loaded onto the Panther system up to 5 times. The Panther system logs each time the reagents are loaded.

Aptima® General Information

D. After thawing the calibrator, the solution must be clear, i.e., not cloudy or have precipitates. Ensure that precipitates are dissolved. Do not use the calibrator if gelling, precipitation, or cloudiness is present.

- - F. The qHBV Target Enhancer Reagent must be at 15°C to 30°C before use.

Specimen Collection and Storage

Note: Handle all specimens as if they contain potentially infectious agents. Use Universal Precautions.

Note: Take care to avoid cross-contamination during sample handling steps. For example, discard used material without passing over open tubes.

Whole blood specimens collected in the following glass or plastic tubes may be used:

- Tubes containing ethylenediaminetetraacetic acid (EDTA) or acid citrate dextrose (ACD) anticoagulants
- · Plasma preparation tubes
- · Serum tubes
- · Serum separator tubes

For serum, allow the clot to form before further processing.

A. Specimen Collection

Whole blood can be stored at 2° C to 30° C for up to 24 hours, and the plasma must be separated via centrifugation in a primary tube prior to processing. Separate the plasma or serum from the red blood cells following the manufacturer's instructions for the tube used. Plasma or serum can be tested on the Panther system in the primary tube or transferred to the secondary Aptima Specimen Aliquot Tube (SAT). To obtain the 500 uL reaction volume, the minimum volume of serum or plasma for primary collection tubes is 1200 μ L and 700 uL for SATs.

If not tested immediately, plasma and serum can be stored in accordance with the specifications below. If transferred to the Specimen Aliquot Tube (SAT), plasma or serum may be frozen at -20°C. Do not exceed 3 freeze—thaw cycles. Do not freeze specimens in EDTA, ACD, or serum primary collection tubes.

B. Specimen Storage Conditions

1. EDTA and ACD Plasma Specimens

Whole blood can be stored at 2°C to 30°C and must be centrifuged within 24 hours of specimen collection. Plasma may then be stored under one of the following conditions:

- In the primary collection tube or SAT at 2°C to 30°C for up to 24 hours,
- In the primary collection tube or SAT at 2°C to 8°C for up to 5 days, or
- In the SAT at -20°C for up to 60 days.

General Information Aptima®

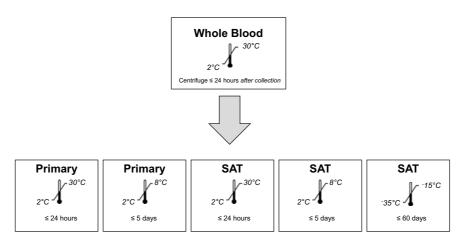


Figure 1. Storage Conditions for EDTA/ACD Tubes

2. PPT Specimens

Whole blood can be stored at 2°C to 30°C and must be centrifuged within 24 hours of specimen collection. Plasma may then be stored under one of the following conditions:

- In the PPT or SAT at 2°C to 30°C for up to 24 hours,
- In the PPT or SAT at 2°C to 8°C for up to 5 days, or
- In the PPT or SAT at -20°C for up to 60 days.

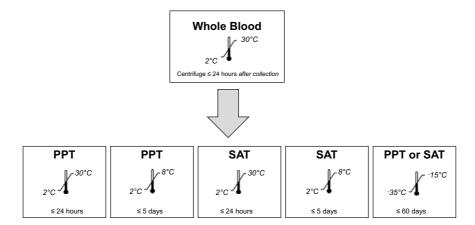


Figure 2. Storage Conditions for PPTs

3. Serum Tube Specimens

Whole blood can be stored at 2°C to 30°C and must be centrifuged within 24 hours of specimen collection. Serum may then be stored under one of the following conditions:

- In the serum tube or SAT at 2°C to 30°C for up to 24 hours,
- In the serum tube or SAT at 2°C to 8°C for up to 5 days, or
- In the SAT at -20°C for up to 60 days.

Aptima® General Information

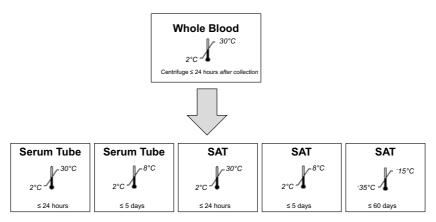


Figure 3. Storage Conditions for Serum Tubes

4. SST Specimens

Whole blood can be stored at 2°C to 30°C and must be centrifuged within 24 hours of specimen collection. Serum may then be stored under one of the following conditions:

- In the SST or SAT at 2°C to 30°C for up to 24 hours,
- In the SST or SAT at 2°C to 8°C for up to 5 days, or
- In the SST or SAT at -20°C for up to 60 days.

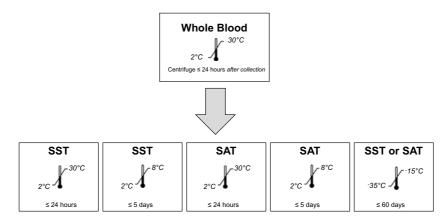


Figure 4. Storage Conditions for SSTs

C. Long Term Frozen Storage

Plasma or serum samples may be stored at -70°C for up to 60 days in SATs.

D. Dilution of Plasma and Serum Specimens

Plasma and serum specimens may be diluted in the SAT for testing on the Panther system. See *Panther System Test Procedure*, step E.7 below for more information.

Note: If a specimen is diluted, it should be tested immediately after dilution. Do not freeze a diluted specimen.

Samples Onboard the Panther System

Samples may be left on the Panther system uncapped for up to 8 hours. Samples may be removed from the Panther system and tested as long as the total time onboard does not exceed 8 hours prior to the pipetting of the sample by the Panther system.

General Information Aptima®

Specimen Transport

Maintain sample storage conditions as described in Specimen Collection and Storage.

Note: Specimens must be shipped in accordance with applicable national, international, and regional transportation regulations.

Aptima® Panther System

Panther System

Reagents for the Aptima HBV Quant assay are listed below for the Panther® system. Reagent identification symbols are also listed next to the reagent name.

Reagents and Materials Provided

Aptima HBV Quant Assay Kit, 100 tests (Cat. No. PRD-03868) (1 assay box, 1 calibrator kit, 1 controls kit, and 1 Target Enhancer Reagent box)

Additional calibrators and controls can be ordered separately. See the individual catalog numbers below.

Aptima HBV Quant Assay Box

(store at 2°C to 8°C upon receipt)

Symbol	Component	Quantity
Α	qHBV Amplification Reagent Non-infectious nucleic acids dried in buffered solution.	1 vial
E	qHBV Enzyme Reagent Reverse transcriptase and RNA polymerase dried in HEPES buffered solution.	1 vial
PRO	qHBV Promoter Reagent Non-infectious nucleic acids dried in buffered solution.	1 vial
AR	qHBV Amplification Reconstitution Solution Aqueous solution containing glycerol and preservatives.	1 x 7.2 mL
ER	qHBV Enzyme Reconstitution Solution HEPES buffered solution containing a surfactant and glycerol.	1 x 5.8 mL
PROR	qHBV Promoter Reconstitution Solution Aqueous solution containing glycerol and preservatives.	1 x 4.5 mL
TCR	qHBV Target Capture Reagent Nucleic acids in a buffered salt solution containing solid phase, non-infectious nucleic acids, and Internal Calibrator.	1 x 72.0 mL
	Reconstitution Collars	3
	Master Lot Barcode Sheet	1 sheet

Aptima HBV Quant Calibrator Kit (Cat. No. PRD-03872)

(store at -15°C to -35°C upon receipt)

Symbol	Component	Quantity
PCAL	qHBV Positive Calibrator Plasmid DNA in buffered solution	5 x 2.5 mL
	Calibrator Barcode Label	_

Aptima HBV Quant Assay 11 AW-15644, Rev. 001

Panther System Aptima®

Aptima HBV Quant Controls Kit (Cat. No. PRD-03869)

(store at -15°C to -35°C upon receipt)

Symbol	Component	Quantity
NC	qHBV Negative Control HBV negative defibrinated human plasma containing gentamicin and 0.2% sodium azide as preservatives.	5 x 0.8 mL
LPC	qHBV Low Positive Control Inactivated HBV positive plasma in defibrinated human plasma containing gentamicin and 0.2% sodium azide as preservatives.	5 x 0.8 mL
НРС	qHBV High Positive Control Inactivated HBV positive plasma in defibrinated human plasma containing gentamicin and 0.2% sodium azide as preservatives.	5 x 0.8 mL
	Control Barcode Label	_

Aptima HBV Quant Target Enhancer Reagent Box

(store at 15°C to 30°C upon receipt)

Symbol	Component	Quantity
TER	qHBV Target Enhancer Reagent A concentrated solution of lithium hydroxide	1 x 46.0 mL

Materials Required but Available Separately

Note: Materials available from Hologic have catalog numbers listed, unless otherwise specified.

Material	Cat. No.
Panther System	_
Aptima HBV Quant Calibrator Kit	PRD-03872
Aptima HBV Quant Controls Kit	PRD-03869
Panther Run Kit for Real Time Assays (for real time assays only)	PRD-03455 (5000 tests)
Aptima Assay Fluids Kit (also known as Universal Fluids Kit) contains Aptima Wash Solution, Aptima Buffer for Deactivation Fluid, and Aptima Oil Reagent	303014 (1000 tests)
Multi-tube units (MTUs)	104772-02
Panther Waste Bag Kit	902731
Panther Waste Bin Cover	504405
Or, Panther System Run Kit	303096 (5000 tests)
(when running non-real time-TMA assays in parallel with real time-TMA assays contains MTUs, waste bags, waste bin covers, auto detect, and assay fluids	s)
Tips, 1000 μL conductive, liquid sensing	10612513 (Tecan)

Aptima® Panther System

Material	Cat. No.
Bleach, 5% to 7% (0.7 M to 1.0 M) sodium hypochlorite solution	_
Disposable, powderless gloves	_
Replacement non-penetrable caps	103036A
Reagent replacement caps Amplification, Enzyme, and Promoter reagent reconstitution bottles TCR bottle TER bottle	CL0041 (100 caps) CL0040 (100 caps) 501604 (100 caps)
Plastic-backed laboratory bench covers	_
Lint-free wipes	_
Pipettor	_
Tips	_
Primary collection tubes of the following dimensions may be used: 13 mm x 100 mm 13 mm x 75 mm 16 mm x 100 mm	_
Centrifuge	_
Vortex mixer	_

Optional Materials

Material	Cat. No.
Aptima Specimen Aliquot Tubes (SAT) (100 pack)	503762
Transport Tube Cap (100 pack) cap for SAT	504415
Aptima Specimen Diluent	PRD-03503
Aptima Specimen Diluent Kit contains specimen diluent, 100 SATs and 100 caps	PRD-03654
Transfer pipets	_
Cotton-tipped swabs	_
Tube rocker	_

Panther System Test Procedure

Note: See the Panther System Operator's Manual for additional procedural information.

A. Work Area Preparation

1. Clean work surfaces where reagents will be prepared. Wipe down work surfaces with 2.5% to 3.5% (0.35 M to 0.5 M) sodium hypochlorite solution. Allow the sodium hypochlorite solution to contact surfaces for at least 1 minute and then follow with a deionized (DI) water rinse. Do not allow the sodium hypochlorite solution to dry. Cover the bench surface with clean, plastic-backed absorbent laboratory bench covers.

Panther System Aptima®

2. Clean a separate work surface where samples will be prepared. Use the procedure described above (step A.1).

3. Clean any pipettors. Use the cleaning procedure described above (step A.1).

B. Calibrator and Controls Preparation

Allow the calibrator and controls to reach 15°C to 30°C prior to processing as follows:

1. Remove the calibrator and controls from storage (-15°C to -35°C) and place at 15°C to 30°C. Throughout the thawing process, gently invert each tube to mix thoroughly. Ensure tube contents are fully thawed prior to use.

Option. Calibrator and control tubes may be placed on a tube rocker to mix thoroughly. Ensure tube contents are fully thawed prior to use.

Note: Avoid creating excessive foam when inverting the calibrator and controls. Foam compromises the level-sensing by the Panther system.

- 2. When the tube contents have thawed, dry the outside of the tube with a clean, dry disposable wipe.
- 3. To prevent contamination, do not open the tubes at this time.
- C. Reagent Reconstitution/Preparation of a New Kit

Note: Reconstitution of reagents should be performed prior to beginning any work on the Panther system.

- 1. To prepare TCR, perform the following:
 - a. Remove the TCR from storage (2°C to 8°C). Check the lot number on the TCR bottle to make sure that it matches the lot number on the Master Lot Barcode Sheet.
 - b. Immediately shake the TCR bottle vigorously 10 times. Allow the TCR bottle to remain at 15°C to 30°C to warm for at least 45 minutes. During this period, swirl and invert the TCR bottle at least every 10 minutes.
 - **Option.** The TCR bottle may be prepared on a tube rocker by following these instructions: Remove the TCR from storage (2°C to 8°C) and immediately shake vigorously 10 times. Place the TCR bottle on a tube rocker and leave the TCR at 15°C to 30°C to warm for at least 45 minutes.
 - c. Ensure all precipitate is in solution and the magnetic particles are suspended before use.
- 2. To reconstitute Amplification, Enzyme, and Promoter Reagents, perform the following:
 - a. Remove the lyophilized reagents and corresponding reconstitution solutions from storage (2°C to 8°C). Pair each reconstitution solution with its lyophilized reagent.
 - b. Ensure that the reconstitution solution and lyophilized reagent have matching label colors. Check the lot numbers on the Master Lot Barcode Sheet to ensure that the appropriate reagents are paired.
 - Open the lyophilized reagent vial by removing the metallic seal and rubber stopper.
 - ii. Firmly insert the notched end of the reconstitution collar (black) onto the vial (Figure 5, Step 1).
 - iii. Open the matching reconstitution solution bottle, and set the cap on a clean, covered work surface.

Aptima® Panther System

iv. Place the reconstitution solution bottle on a stable surface (i.e., bench). Then, invert the lyophilized reagent vial over the reconstitution solution bottle and firmly attach the collar to the reconstitution solution bottle (Figure 5, Step 2).

- v. Slowly invert the assembled bottles (vial attached to solution bottle) to allow the solution to drain into the glass vial (Figure 5, Step 3).
- vi. Pick up the assembled bottles, and swirl the assembled bottles for at least 10 seconds (Figure 5, Step 4).
- vii. Wait at least 30 minutes for the lyophilized reagent to go into solution.
- viii. After the lyophilized reagent has gone into solution, swirl the assembled bottles for at least 10 seconds and then slightly rock the solution within the glass vial back and forth to mix thoroughly.
- c. Slowly tilt the assembled bottles again to allow all of the solution to drain back into the reconstitution solution bottle (Figure 5, Step 5).
- d. Carefully remove the reconstitution collar and glass vial (Figure 5, Step 6).
- e. Recap the bottle. Record operator initials and reconstitution date on the label (Figure 5, Step 7).
- f. Discard the reconstitution collar and glass vial (Figure 5, Step 8).

Warning: Avoid creating excessive foam when reconstituting reagents. Foam compromises the level-sensing by the Panther system.

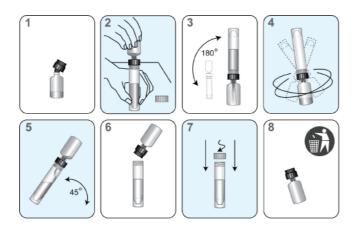


Figure 5. Reagent Reconstitution Process

- 3. Remove the qHBV Target Enhancer Reagent from storage (15°C to 30°C). Record operator initials and open date on the label. Check the lot number on the TER bottle to make sure it matches the lot number on the Master Lot Barcode Sheet.
- D. Reagent Preparation for Previously Prepared Reagents
 - 1. Remove the previously prepared reagents from storage (2°C to 8°C). Previously prepared Amplification, Enzyme and Promoter reagents, and TCR must reach 15°C to 30°C prior to the start of the assay.
 - 2. Remove TER from storage (15°C to 30°C).
 - 3. For previously prepared TCR, perform step C.1 above prior to loading on the system.
 - 4. Swirl and invert the Amplification, Enzyme, and Promoter reagents to mix thoroughly prior to loading on the system. Avoid creating excessive foam when inverting reagents.

Panther System Aptima®

5. Do not top off reagent bottles. The Panther system will recognize and reject bottles that have been topped off.

E. Specimen Handling

- Ensure that processed specimens in primary tubes or undiluted plasma specimens in SATs have been stored properly per Specimen Collection and Storage section. Additional centrifugation is only recommended at any time when plasma and red blood cells have mixed again in the primary tube, or when bubbles are present. Bubbles can compromise the level sensing by the Panther system.
- 2. Ensure frozen specimens are thoroughly thawed.
- 3. Allow the specimens to reach 15°C to 30°C prior to processing. See Samples Onboard the Panther System for additional onboard information.
- 4. Ensure each primary collection tube contains at least 1200 μL of specimen. Ensure each Aptima SAT contains at least 700 μL of specimen. If specimen dilution is necessary, see step E.7 below for additional information.
- 5. Vortex previously frozen specimens in SATs for 3 to 5 seconds to mix thoroughly and prior to loading into a sample rack, centrifuge at 1000 to 3000g for 10 minutes.
- 6. If specimen centrifugation is required, centrifuge each specimen at 1000 to 3000g for 10 minutes just prior to loading specimens into a Sample Rack. Do not remove caps. Bubbles in the tube compromise the level-sensing by the Panther system.
 - See *System Preparation*, step F.2 below, for information about loading the rack and removing the caps.
- 7. Diluting a plasma or serum specimen in the SAT:

A specimen may be diluted in the SAT for testing on the Panther system.

Note: If a specimen is diluted, it must be tested immediately after dilution.

a. Dilution of low-volume specimens

The volume of specimens may be increased to the minimum volume required (700 μ L) using Aptima Specimen Diluent. Specimens with at least 240 μ L may be diluted with two parts specimen diluent (1:3) as follows:

- i. Place 240 µL of specimen in the SAT.
- ii. Add 480 µL of Specimen Diluent.
- iii. Cap the tube.
- iv. Gently invert 5 times to mix.

Specimens diluted 1:3 can be tested using the 1:3 option on the Panther system (see the *Panther System Operator's Manual* for more information). The software will automatically report the neat result by applying the dilution factor. These specimens will be flagged as diluted specimens.

b. Dilution of high-titer specimens

If a specimen's result is above the upper limit of quantitation (ULoQ), it may be diluted with 99 parts of Aptima Specimen Diluent (1:100) as follows:

- i. Place 30 µL of specimen in the SAT.
- ii. Add 2970 µL of specimen diluent.
- iii. Cap the tube.
- iv. Gently invert 5 times to mix.

Specimens diluted 1:100 can be tested using the 1:100 option on the Panther system (see *Panther System Operator's Manual* for more information). The

Aptima® Panther System

software will automatically report the neat result by applying the dilution factor. These specimens will be flagged as diluted specimens.

Note: For diluted specimens with neat concentrations greater than the ULoQ, results will be reported using scientific notation.

F. System Preparation

- 1. Set up the system according to the instructions in the *Panther System Operator's Manual* and *Procedural Notes*. Make sure that the appropriately sized reagent racks and TCR adapters are used.
- 2. Load samples into the Sample Rack. Perform the following steps for each sample tube (specimen and, when necessary, calibrator and controls):
 - a. Loosen one sample tube cap, but do not remove it yet.

Note: Be especially careful to avoid contamination by the spread of aerosols. Gently loosen caps on samples.

- b. Load the sample tube into the Sample Rack.
- c. Repeat steps 2.a and 2.b for each remaining sample.
- d. After the samples have been loaded into the Sample Rack, remove and discard each sample tube cap in one Sample Rack. To avoid contamination, do not pass a cap over any other Sample Racks or sample tubes.
- e. If necessary, use a new, disposable transfer pipet to remove any bubbles or foam.
- f. When the last cap has been removed, load the Sample Rack into the Sample Bay.

Note: If running other assays and sample types at the same time, secure the Sample Retainer prior to loading the Sample Rack into the Sample Bay.

g. Repeat steps 2.a to 2.f for the next Sample Rack.

Procedural Notes

A. Calibrator and Controls

- 1. The qHBV positive calibrator, qHBV low positive control, qHBV high positive control, and qHBV negative control tubes can be loaded in any position in the Sample Rack and in any Sample Bay Lane on the Panther system. Specimen pipetting will begin when one of the following two conditions has been met:
 - a. The calibrator and controls are currently being processed by the system.
 - b. Valid results for the calibrator and controls are registered on the system.
- Once the calibrator and control tubes have been pipetted and are processing for the Aptima HBV Quant assay reagent kit, specimens can be tested with the associated reconstituted kit for up to 24 hours unless:
 - a. The calibrator or control results are invalid.
 - b. The associated assay reagent kit is removed from the system.
 - c. The associated assay reagent kit has exceeded stability limits.
- 3. The calibrator and each control tube can be used once. Attempts to use the tube more than once can lead to processing errors.

B. Glove Powder

As in any reagent system, excess powder on some gloves may cause contamination of opened tubes. Powderless gloves are recommended.

Quality Control Aptima®

Quality Control

A run or specimen result may be invalidated by an operator if technical, operator, or instrument difficulties are observed while performing the assay and are documented. In this case, specimens must be retested.

Assay Calibration

To generate valid results, an assay calibration must be completed. A single positive calibrator is run in triplicate each time a reagent kit is loaded on the Panther system. Once established, the calibration is valid for up to 24 hours. Software on the Panther system alerts the operator when a calibration is required. The operator scans a calibration coefficient found on the Master Lot Barcode Sheet provided with each reagent kit.

During processing, criteria for acceptance of the calibrator are automatically verified by the software on the Panther system. If less than two of the calibrator replicates are valid, the software automatically invalidates the run. Samples in an invalidated run must be retested using a freshly prepared calibrator and freshly prepared controls.

Negative and Positive Controls

To generate valid results, a set of assay controls must be tested. One replicate of the negative control, the low positive control, and the high positive control must be tested each time a reagent kit is loaded on the Panther system. Once established, the controls are valid for up to 24 hours. Software on the Panther system alerts the operator when controls are required.

During processing, criteria for acceptance of controls are automatically verified by software on the Panther system. To generate valid results, the negative control must give a result of "Not Detected" and the positive controls must give results within predefined parameters (LPC Nominal Target: 2.7 log₁₀ IU/mL, HPC Nominal Target: 4.6 log₁₀ IU/mL. If any one of the controls has an invalid result, the software automatically invalidates the run. Samples in an invalidated run must be retested using a freshly prepared calibrator and freshly prepared controls.

Internal Calibrator/Internal Control

Each sample contains an internal calibrator/internal control (IC). During processing, IC acceptance criteria are automatically verified by the Panther system software. If an IC result is invalid, the sample result is invalidated. Every sample with an invalid IC result must be retested to obtain a valid result.

The Panther system software is designed to accurately verify processes when procedures are performed following the instructions provided in this package insert and the *Panther System Operator's Manual.*

Interpretation of Results

The Panther system automatically determines the concentration of HBV DNA for specimens and controls by comparing the results to a calibration curve. HBV DNA concentrations are reported in IU/mL and log₁₀ IU/mL. The interpretation of results is provided in Table 1. If the dilution option is used, the Panther system automatically calculates the HBV DNA concentration for the neat specimen by multiplying the diluted concentration by the dilution factor, and diluted samples will be flagged as diluted.

Note: For diluted specimens, results listed as "Not Detected" or "<10 detected" may be generated by diluting a specimen with a concentration above, but close to the LoD (limit of detection) or LLoQ (lower limit of quantitation). It is recommended to collect and test another neat specimen if a quantitative result is not obtained

Table 1: Result Interpretation

Reported Aptima HBV Quant Assay Result		l
IU/mL	Log ₁₀ IU/mL ^a	- Interpretation
Not Detected	Not Detected	HBV DNA not detected.
<10 detected	<1.00	HBV DNA is detected but at a level below the LLoQ
10 to 1,000,000,000	1.00 to 9.00	HBV DNA concentration is within the linear range of 10 to 1,000,000,000 IU/mL
> 1,000,000,000	> 9.00	HBV DNA concentration is above the ULoQ
Invalid ^b	Invalid ^b	Error indicated in the generation of the result. Specimen should be retested

^a Value is truncated to two decimal places.

Note: For diluted specimens with neat concentrations greater than the ULoQ, results will be reported using scientific notation.

Limitations

- A. Use of this assay is limited to personnel who have been trained in the procedure. Failure to follow the instructions given in this package insert may result in erroneous results.
- B. Reliable results are dependent on adequate specimen collection, transport, storage, and processing.

^b Invalid results are displayed in blue colored font.

Nonclinical Performance

Limit of Detection Using the 3rd WHO International Standard

The limit of detection (LoD) of the assay is defined as the concentration of HBV DNA that is detected at 95% or greater probability according to CLSI EP17-A2.¹²

The LoD was determined by testing panels of the 3rd WHO International Standard for Hepatitis B Virus DNA (NIBSC 10/264, genotype A) diluted in HBV negative human EDTA plasma and serum. A minimum of 36 replicates of each dilution were tested with each of three reagent lots for a minimum of 108 replicates per dilution. Probit analysis was performed to generate the 95% predicted detection limits. The LoD values are the results from the reagent lot with the highest predicted detection limit. The LoD for the Aptima HBV Quant assay using the 3rd WHO International Standard is 4.8 IU/mL for plasma and 5.9 IU/mL for serum.

Limit of Detection Across HBV Genotypes

The LoD was determined by testing dilutions of HBV positive clinical specimens for genotypes A, B, C, D, E, F, G, and H in HBV negative human plasma and serum. Concentrations were determined using an FDA approved assay. A minimum of 24 replicates of each panel member were tested with each of two reagent lots for a minimum of 48 replicates per panel member. Probit analysis was performed to generate 95% predicted detection limits. The LoD values shown in Table 2 are the results from the reagent lot with the highest predicted detection limit.

Table 2: Limit of Detection (95% predicted detection limit) Across HBV Genotypes Using Clinical Specimens

_	Concentration (IU/mL)		
Genotype	Plasma	Serum	
А	3.3	4.1	
В	2.9	3.9	
С	4.9	5.2	
D	5.7	5.4	
Е	5.8	5.8	
F	3.0	4.0	
G	2.8	7.4	
Н	5.5	6.3	

Linear Range

The linear range was established by testing panels of HBV genotype A virus (0.78 \log_{10} IU/mL to 7.30 \log_{10} IU/mL) and plasmid DNA (5.30 \log_{10} IU/mL to 9.18 \log_{10} IU/mL) diluted in HBV negative human plasma and serum according to CLSI EP06-A.¹³ The Aptima HBV Quant assay demonstrated linearity across the range tested, with an upper limit of quantitation (ULoQ) of 9 \log_{10} IU/mL as shown in Figure 6.

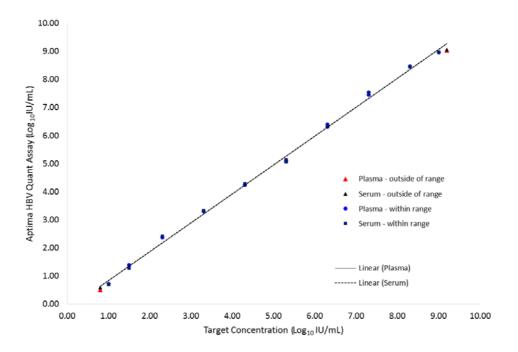


Figure 6. Linearity in Plasma and Serum

Linearity Across HBV Genotypes

The linearity of HBV genotypes was established by testing individual clinical positive samples for genotypes A, E, F, G, and H, and PEI 1st WHO Reference panels (PEI 5086/08) for genotypes B, C, and D. Virus was used for the lower range of the assay (4 \log_{10} IU/mL and below for genotypes A-G, 3 \log_{10} IU/mL, and below for genotype H). Plasmid DNA was used for the upper range with a 2 \log_{10} overlap. Dilutions in negative human plasma were tested for all genotypes. Linearity was demonstrated across the range tested for all genotypes tested as shown in Figure 7.

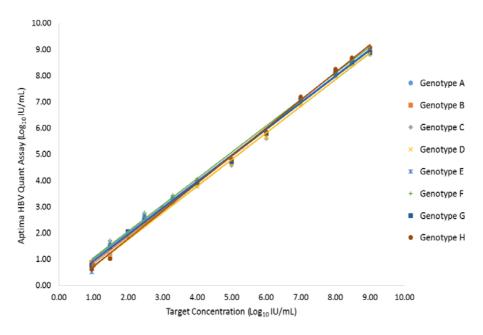


Figure 7. Linear Range and Linearity (Plasma)

Lower Limit of Quantitation Using the 3rd WHO International Standard

The lower limit of quantitation is defined as the lowest concentration at which HBV DNA is reliably quantitated within a total error, according to CLSI EP17-A2.¹² Total error was estimated by the following two methods:

Total Analytical Error (TAE) = |bias| + 2SD and Total Error (TE) = SQRT(2) x 2SD.

To ensure accuracy and precision of measurements, the total error of the Aptima HBV Quant assay was set at $1 \log_{10} IU/mL$ (i.e., at the LLoQ, the difference between two measurements of more than $1 \log_{10} IU/mL$ is statistically significant).

The LLoQ was determined by testing panels of the 3rd WHO International Standard for Hepatitis B Virus DNA (NIBSC 10/264, genotype A) diluted in HBV negative human plasma and serum. Forty-five (45) replicates of each dilution were tested with each of three reagent lots for a minimum of 135 replicates per dilution. The results for the three reagent lots are shown in Table 3 for plasma and Table 4 for serum. The results for the lowest observed concentration that met the accuracy goal (TE \leq 1 log₁₀ IU/mL and TAE \leq 1 log₁₀ IU/mL) with >95% detection and greater than or equal to the LoD are shaded in both tables and summarized in Table 5.

The calculated LLoQ for the 3rd WHO International Standard for Hepatitis B Virus is 6 IU/mL (0.79 log₁₀IU/mL) for plasma and 8 IU/mL (0.88 log₁₀IU/mL) for serum, which are based on the highest calculated concentration among the three reagent lots in accordance with EP 17-A2. The LLoQ was established across genotypes (see next section "Lower Limit of Quantitation Across HBV Genotypes"). This genotype data establishes the overall LLoQ for the assay as 10 IU/mL.

Table 3: Determination of LLoQ Using the 3rd WHO International Standard for HBV Diluted in Plasma

Reagent Lot	Aptima HBV Quant (IU/mL)	Aptima HBV Quant (Log ₁₀ IU/mL)	SD (Log ₁₀ IU/mL)	Bias (Log ₁₀ IU/mL)	Calculated TE (Log ₁₀ IU/mL)	Calculated TAE (Log ₁₀ IU/mL)
	3	0.53	0.21	0.32	0.59	0.74
1	3	0.52	0.21	0.38	0.61	0.81
Ī	5	0.70	0.23	0.25	0.65	0.71
	6	0.76	0.26	0.09	0.73	0.60
2	5	0.69	0.22	0.21	0.63	0.65
	6	0.77	0.25	0.18	0.70	0.68
3	6	0.79	0.31	0.05	0.88	0.68
	8	0.88	0.23	0.02	0.66	0.48
	9	0.96	0.23	0.00	0.66	0.47

SD=standard deviation

Panel members that met the accuracy goal (TE ≤ 1 and TAE ≤ 1), > LoD, and > 95% detection for Reagent Lots 1, 2, and 3 are shaded.

Table 4: Determination of LLoQ Using the 3rd WHO International Standard for HBV Diluted in Serum

Reagent Lot	Aptima HBV Quant (IU/mL)	Aptima HBV Quant (Log ₁₀ IU/mL)	SD (Log ₁₀ IU/mL)	Bias (Log ₁₀ IU/mL)	Calculated TE (Log ₁₀ IU/mL)	Calculated TAE (Log ₁₀ IU/mL)
	4	0.65	0.24	0.20	0.67	0.67
1	5	0.65	0.24	0.25	0.68	0.73
_	5	0.67	0.22	0.28	0.63	0.73
-	5	0.70	0.29	0.14	0.82	0.72
2	5	0.72	0.27	0.18	0.77	0.72
_	6	0.75	0.24	0.20	0.68	0.68
	8	0.88	0.29	0.04	0.83	0.63
3	10	0.98	0.23	0.08	0.66	0.55
_	11	1.02	0.28	0.07	0.78	0.62

SD=standard deviation

Panel members that met the accuracy goal (TE ≤ 1 and TAE ≤ 1), > LoD, and > 95% detection for Reagent Lots 1, 2, and 3 are shaded.

Table 5: Summary of the Calculated LLoQ Using the 3rd WHO International Standard for HBV

Dogwood Lot	Plas	ma LLoQ	Se	rum LLoQ
Reagent Lot	(IU/mL)	(Log ₁₀ IU/mL)	(IU/mL)	(Log ₁₀ IU/mL)
1	5	0.70	4	0.65
2	5	0.69	5	0.72
3	6	0.79	8	0.88

Lower Limit of Quantitation Across HBV Genotypes

The LLoQ was determined by testing dilutions of HBV positive clinical specimens for genotypes A, B, C, D, E, F, G, and H in HBV negative human plasma and serum. Assignment of the concentration for clinical specimens was determined using an FDA-approved assay. A total of 36 replicates of each panel member were tested with each of two reagent lots for a minimum of 72 replicates per panel member. The results for the lowest observed concentration meeting the accuracy goal (TE \leq 1 log₁₀ IU/mL and TAE \leq 1 log₁₀ IU/mL) with >95% detection for each reagent lot are shown in Table 6 for plasma and Table 7 for serum. The highest observed concentration among the reagent lots for each genotype is summarized in Table 8. Genotype D in serum had the highest LLoQ at 9 IU/mL (0.96 log₁₀ IU/mL). This supported the overall LLoQ for the assay as 10 IU/mL.

Table 6: Determination of LLoQ Across HBV Genotypes in Plasma

HBV Genotype	Reagent Lot	Aptima HBV Quant ^a (IU/mI)	Aptima HBV Quant ^a (Log ₁₀ IU/ml)	SD (Log ₁₀ IU/ml)	Bias (Log ₁₀ IU/ml)	Calculated TE (Log ₁₀ IU/ml)	Calculated TAE (Log ₁₀ IU/ml)
A	1	6	0.74	0.24	0.10	0.67	0.58
A	2	7	0.85	0.20	0.00	0.57	0.40
В	1	4	0.65	0.24	0.20	0.67	0.67
Ь	2	6	0.75	0.25	0.10	0.70	0.59
С	1	4	0.61	0.21	0.35	0.60	0.77
C	2	6	0.75	0.30	0.20	0.84	0.80
	1	7	0.87	0.28	0.31	0.81	0.88
D	2	8	0.91	0.30	0.26	0.86	0.87
	1	8	0.88	0.32	0.29	0.89	0.92
E	2	7	0.82	0.21	0.36	0.60	0.78
	1	6	0.76	0.27	0.08	0.76	0.62
F	2	7	0.86	0.30	0.01	0.84	0.60
	1	4	0.59	0.21	0.25	0.60	0.68
G	2	4	0.65	0.23	0.19	0.64	0.64
	1	6	0.78	0.28	0.39	0.80	0.96
Н	2	7	0.83	0.27	0.34	0.78	0.89

SD=standard deviation

^a Additional levels were run, but that data is not reported in the table.

Table 7: Determination of LLoQ Across HBV Genotypes in Serum

HBV Genotype	Reagent Lot	Aptima HBV Quant ^a	Aptima HBV Quant ^a	SD	Bias	Calculated TE	Calculated TAE
,	Ū	(IU/ml)	(Log ₁₀ IU/ml)	(Log ₁₀ IU/ml)	(Log ₁₀ IU/ml)	(Log ₁₀ IU/ml)	(Log ₁₀ IU/ml)
A	1	4	0.65	0.26	0.25	0.73	0.77
A	2	6	0.81	0.26	0.04	0.74	0.56
	1	5	0.70	0.19	0.26	0.54	0.64
В	2	5	0.72	0.25	0.18	0.72	0.69
-	1	5	0.66	0.26	0.30	0.74	0.82
С	2	6	0.81	0.26	0.10	0.74	0.62
	1	7	0.84	0.27	0.42	0.75	0.95
D	2	9	0.96	0.27	0.29	0.76	0.83
	1	6	0.79	0.26	0.38	0.75	0.91
E	2	8	0.89	0.28	0.29	0.79	0.84
	1	6	0.76	0.23	0.08	0.66	0.55
F	2	5	0.71	0.26	0.14	0.72	0.65
	1	8	0.89	0.28	0.29	0.80	0.85
G	2	5	0.66	0.24	0.29	0.67	0.77
	1	6	0.74	0.24	0.43	0.67	0.90
Н	2	6	0.81	0.29	0.37	0.82	0.95

SD=standard deviation

Table 8: Summary of LLoQ Across Genotypes in Plasma and Serum

Ganatyna	Plas	ma LLoQ	Serum LLoQ		
Genotype	(IU/mL) (Log ₁₀ IU/mL)		(IU/mL)	(Log ₁₀ IU/mL)	
Α	7	0.85	6	0.81	
В	6	0.75	5	0.72	
С	6	0.75	6	0.81	
D	8	0.91	9	0.96	
E	8	0.88	8	0.89	
F	7	0.86	6	0.76	
G	4	0.65	8	0.89	
Н	7	0.83	6	0.81	

^a Additional levels were run, but that data is not reported in the table.

Traceability to the 3rd WHO International Standard

A series of secondary standards with known concentrations was used throughout product development and product manufacturing to establish traceability to the WHO standard. The concentrations tested for the HBV WHO standard were between 2.0 and 4.0 log₁₀ IU/mL, the secondary standards ranged in concentration from 2.4 to 8.4 log₁₀ IU/mL. The Aptima HBV Quant assay controls and calibrators were also tested along with the secondary standards and the WHO standard. All of the panels had similar results, and they were distributed linearly across the assay's linear range, as presented in Figure 8.

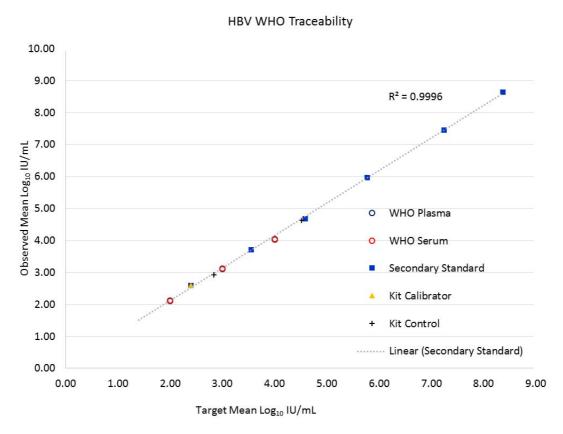


Figure 8. Traceability Between the 3rd HBV WHO Standard Target Concentrations and Observed Concentration From the Aptima HBV Quant Assay

Precision

The Aptima HBV Quant precision panel was built by diluting HBV genotype A virus and HBV plasmid DNA into HBV-negative clinical plasma and HBV-negative clinical serum (the four highest panel members in each matrix were plasmid DNA). Eleven panel members in each matrix spanned the range of the assay (target concentrations of 1.30 log₁₀ IU/mL to 8.90 log₁₀ IU/mL) and were tested in three replicates per run by one operator, using three reagent lots on one Panther system over three days, two runs a day.

Table 9 shows the precision of assay results (in log_{10} IU/mL) between days, between lots, between runs, within runs, and overall. Total variability was primarily due to the intra-run measurement (i.e., random error).

Aptima® **Nonclinical Performance**

Table 9: Precision of the Aptima HBV Quant Assay

NA4	0		Mean Concentration	Mean Concentration	Inter-Lot	Inter-Lot	Inter-Day	Inter-Day	Inter-Run	Inter-Run	Intra-Run	Intra-Run	Total	Total
Matrix	Sample	N	(IU/mL)	(log ₁₀ IU/mL)	SD	Log Normal CV (%)	SD	Log Normal CV (%)						
		34 ^a	32	1.21	0.07	16.2	0.07	16.2	0.05	11.7	0.28	71.8	0.3	78.2
		54	97	1.95	0.05	11.7	0.03	6.8	0.02	5.7	0.15	36.8	0.17	39.9
	Virus	54	1,474	3.16	0	1.0	0.02	3.8	0.02	4.8	0.07	15.6	0.07	16.8
		54	10,602	4.02	0.02	4.8	0.02	3.6	0.01	2.4	0.06	14.9	0.07	16.2
		54	429,428	5.63	0.03	6.7	0.01	2.7	0.01	2.4	0.06	14.1	0.07	16.1
Plasma	Plasmid DNA	54	652,103	5.8	0.06	14.2	0.01	3.4	0.01	2.0	0.06	13.1	0.09	19.8
	Virus	54	7,617,612	6.88	0.02	5.6	0.00	0.4	0.02	4.2	0.06	13.4	0.07	15.1
	Plasmid DNA	54	10,662,942	7.02	0.02	5.3	0.01	2.3	0.01	2.6	0.06	13.1	0.06	14.5
	Virus	54	89,149,358	7.95	0.01	3.4	0.01	2.7	0.01	1.6	0.04	9.7	0.05	10.8
	Plasmid DNA	54	103,400,000	8.01	0.04	9.4	0.02	3.7	0.01	2.1	0.05	12.0	0.07	15.9
	Plasifilia DINA	54	612,200,000	8.78	0.03	7.6	0.01	1.9	0.01	1.5	0.04	9.0	0.05	12.0
		29 ^a	33	1.27	0.13	31.0	0.08	18.6	0.06	13.9	0.29	75.0	0.33	88.4
		54	88	1.92	0.05	11.2	0.02	4.8	0.02	5.5	0.12	29.1	0.14	32.3
	Virus	54	1,446	3.15	0.02	3.7	0.01	1.9	0.00	1.1	0.08	17.8	0.08	18.3
		54	7,873	3.89	0.02	4.8	0.01	2.6	0.01	2.1	0.06	13.4	0.06	14.6
		54	313,518	5.49	0.01	3.3	0.01	3.0	0.01	3.2	0.08	17.4	0.08	18.3
Serum	Plasmid DNA	54	599,225	5.77	0.04	8.5	0.01	2.5	0.01	2.5	0.06	14.7	0.08	17.4
	Virus	54	7,011,440	6.84	0.02	3.5	0.01	3.2	0.01	3.3	0.07	16.9	80.0	17.9
	Plasmid DNA	54	8,845,332	6.94	0.05	11.9	0.01	2.2	0.01	2.2	0.05	12.2	0.07	17.4
	Virus	54	70,350,774	7.84	0.03	6.3	0.02	3.5	0.02	4.4	0.06	14.4	0.07	16.7
	Plasmid DNA	53 ^b	122,800,000	8.08	0.04	9.9	0.01	2.2	0.01	1.8	0.04	10.4	0.06	14.6
	FIASIIIIU DINA	54	678,700,000	8.83	0.02	3.6	0.01	2.2	0.01	1.8	0.05	11.7	0.05	12.5

SD=standard deviation

^bOne replicate had an invalid result. Log Normal CV(%) = sqrt(10^(SD^2 * ln(10)) - 1) * 100

^aDetected replicates with quantifiable result, total detected replicates = 54.

Potentially Interfering Substances

The susceptibility of the Aptima HBV Quant assay to interference by elevated levels of endogenous substances or by drugs commonly prescribed to HBV infected individuals was evaluated. HBV negative plasma samples and samples spiked with HBV to concentrations of approximately 30 IU/mL (1.48 log₁₀ IU/mL) and 20,000 IU/mL (4.30 log₁₀ IU/mL) were tested.

No interference in the performance of the assay was observed in the presence of albumin (90 mg/mL), hemoglobin (5 mg/mL), triglycerides (30 mg/mL), or unconjugated bilirubin (0.2 mg/mL).

Clinical plasma specimens from patients with elevated levels of defined substances or from patients with the diseases, ten samples for each substance, listed in Table 10 were tested with the Aptima HBV Quant assay. No interference in the performance of the assay was observed.

Table 10: Tested Clinical Specimen Types

	Clinical Specimen Types
1	Antinuclear antibody (ANA)
2	Rheumatoid factor (RF)
3	Alcoholic cirrhosis (AC)
4	Alcoholic hepatitis
5	Non-alcoholic hepatitis
6	Autoimmune hepatitis
7	Elevated alanine aminotransferase (ALT)
8	Hepatocellular carcinoma (HCC)
9	Multiple sclerosis (MS)
10	Systemic lupus erythematosus (SLE)
11	Hyperglobulinemia
12	Rheumatoid arthritis (RA)
13	Anti-Jo1 antibody (JO-1)
14	Multiple myeloma (MM)
15	Hemolyzed (elevated hemoglobin)
16	Icteric (elevated bilirubin)
17	Lipemic (elevated lipid)
18	Elevated protein

No interference in the performance of the assay was observed in the presence of the exogenous substances listed in Table 11 at concentrations at least three times the C_{max} (human plasma).

Table 11: Exogenous Substances

Exogenous Substance Pool	Exogenous Substances Tested
1	Saquinavir, ritonavir, amprenavir, indinavir, lopinavir, nelfinavir mesylate
2	Clarithromycin, valganciclovir hydrochloride, efavirenz, nevirapine
3	Paroxetine HCI, enfuvirtide, zidovudine, didanosine, abacavir sulfate
4	Ribavirin, entecavir, adefovir dipivoxil, tenofovir disoproxil fumarate, lamivudine, ganciclovir, acyclovir
5	Stavudine, ciprofloxacin, fluoxetine, azithromycin, valacyclovir, sertraline, zalcitabine
6	Interferon alpha -2a, interferon alpha -2b, pegylated interferon alpha-2b

Analytical Specificity

Potential cross-reactivity to the pathogens listed in Table 12 was evaluated in HBV negative human plasma in the presence or absence of 30 IU/mL (1.5 \log_{10} IU/mL) and 20,000 (4.3 \log_{10} IU/mL) HBV DNA. No cross-reactivity was observed. No interference was observed in the presence of the pathogens.

Table 12: Pathogens Tested for Analytical Specificity

Microorganism/Pathogen	Con	centration	Microorganism/Pathogen	Concentration		
Adenovirus 5	100,000	TCID50/mLa	Candida albicans	1,000,000	CFU/mL	
BK human polyomavirus	1,000	TCID50/mL	Chlamydia trachomatis	1,000,000	IFU/mL ¹	
Cytomegalovirus	100,000	TCID50/mL	Corynebacterium diphtheriae	1,000,000	CFU/mL	
Dengue virus 1	10,000	TCID50/mL	Neisseria gonorrhoeae	1,000,000	CFU/mL	
Dengue virus 2	10,000	TCID50/mL	Propionibacterium acnes	1,000,000	CFU/mL	
Dengue virus 3	10,000	TCID50/mL	Staphylococcus aureus	1,000,000	CFU/ml	
Dengue virus 4	100,000	TCID50/mL	Staphylococcus epidermidis	1,000,000	CFU/mL	
Epstein-Barr virus	100,000	copies/mL	Streptococcus pneumoniae	1,000,000	CFU/mL	
Flu H1N1	100,000	TCID50/mL	Trichomonas vaginalis	1,000,000	CFU/mL	
Hepatitis A virus	100,000	TCID50/mL	a TCID50/mL = Tissue culture Infectious of	lose units per	mL	
Hepatitis C virus	100,000	IU/mL ^b	 b IU/mL = International units per mL			
Hepatitis G virus	100,000	copies/mL	c vp/mL = Viral particles per mL			
Human herpes virus 6B	100,000	copies/mL	d LD50/mL = Lethal dose per mL			
Human herpes virus 8	100,000	copies/mL	 e CFU/mL = Colony forming units per mL			
HIV-1	100,000	IU/mL				
HIV-2	10,000	TCID50/mL	=			
Human papillomavirus	100,000	copies/mL	_			
Herpes simplex virus 1 (HSV-1)	100,000	TCID50/mL	_			
Herpes simplex virus 2 (HSV-2)	100,000	TCID50/mL	_			
Human T-cell lymphotropic virus 1 (HTLV-1)	100,000	vp/mL ^c	_			
Human T-cell lymphotropic virus 2 (HTLV-2)	100,000	vp/mL	_			
Japanese encephalitis virus	N/A	N/A	_			
Murray Valley encephalitis virus	2,000	LD50/mL ^d	_			
Parvovirus B19	100,000	IU/mL	_			
Rubella virus	10,000	TCID50/mL	-			
St. Louis encephalitis virus	100,000	TCID50/mL	-			
Vaccinia virus	1,000	TCID50/mL	_			
West Nile virus	100,000	TCID50/mL	_			
Yellow fever virus	100,000	TCID50/mL	_			

Matrix Equivalency

One hundred eighteen sample sets of matched blood collection tubes (serum tube, ACD, K2 EDTA, K3 EDTA, PPT, SST) were assessed for matrix equivalency. Of these, 44 sets were naturally infected HBV-positive, and 74 sets were HBV-negative spiked with HBV virus. Correlation for each blood collection tube type, as measured using the serum collection tube as a comparator, is shown in Table 14.

Table 13: Matrix Equivalency Study

		95% CI	of Slope	95% CI of	Intercept		
Blood Collection Tube	Deming Regression	Lower Limit	Upper Limit	Lower Limit	Upper Limit	R ²	Mean Difference (Log ₁₀)
ACD	y = 1.01x - 0.04	1.00	1.02	-0.10	0.01	0.998	-0.01
K2 EDTA	y = 1.02x - 0.14	1.00	1.03	-0.20	-0.07	0.997	-0.07
K3 EDTA	y = 1.01x - 0.12	1.00	1.03	-0.18	-0.06	0.997	-0.06
PPT	y = 1.02x - 0.14	1.00	1.03	-0.21	-0.07	0.996	-0.06
SST	y = 1.00x - 0.03	0.99	1.01	-0.07	0.03	0.999	-0.01

CI=confidence interval

Sample Dilution Using Aptima Specimen Diluent (1:3)

To assess detection accuracy of HBV DNA in samples diluted with Aptima Specimen Diluent, samples that spanned the linear range were diluted 1:3 with Aptima Specimen Diluent (such as 240 μ L of sample combined with 480 μ L of Aptima Specimen Diluent). Each sample was tested neat and diluted (1:3) in triplicate. Testing was performed using one lot of assay reagents on two Panther systems with two Aptima Specimen Diluent lots. The differences between the average reported concentration in native matrix (dilution factor applied to the diluted sample result) and the average concentration in Aptima Specimen Diluent are shown in Table 14 for plasma and Table 15 for serum. The sample concentrations were accurately recovered in the diluted samples.

Table 14: Plasma Specimen 1:3 Dilution Matrix Comparison Summary

Plasma Matrix Average Reported Concentration (Log ₁₀ IU/mL) n = 9	Diluent Average Reported Concentration (Log ₁₀ IU/mL) n = 18	Difference Diluent from Plasma Matrix (Log ₁₀ IU/mL)		
1.20 ^a	1.11 ^b	-0.09		
1.56 ^a	1.36 ^b	-0.20		
2.15	2.04	-0.11		
3.10	2.97	-0.13		
3.92	3.89	-0.03		
4.82	4.79	-0.03		
5.70	5.70	0.00		
7.07	6.98	-0.09		
7.74	7.60	-0.14		
8.74	8.62	-0.12		
9.29	9.19	-0.10		
9.39	9.29	-0.10		

^an=21

^bn=42

Table 15: Serum Specimen 1:3 Dilution Matrix Comparison Summary

Serum Matrix Average Reported Concentration (Log ₁₀ IU/mL) n = 9	Diluent Average Reported Concentration (Log ₁₀ IU/mL) n = 18	Difference Diluent from Serum Matrix (Log ₁₀ IU/mL)		
1.21 ^a	1.11 ^b	-0.10		
1.54 ^a	1.36 ^b	-0.18		
2.21	2.03	-0.18		
3.06	2.98	-0.08		
3.90	3.83	-0.07		
4.77	4.76	-0.01		
5.77	5.74	-0.03		
7.03	7.00	-0.03		
7.85	7.71	-0.14		
8.87	8.76	-0.11		
9.37	9.30	-0.07		
9.46	9.36	-0.10		

^an=21

Sample Dilution Using Aptima Specimen Diluent (1:100)

To assess detection accuracy of HBV DNA in samples diluted with Aptima Specimen Diluent, plasma or serum, eight individual plasma specimens and eight individual serum specimens spiked with HBV virus targeting between 6 to 8 log₁₀ IU/mL, along with eight individual plasma specimens and eight individual serum specimens spiked with HBV plasmid DNA targeting 9.16 log₁₀ IU/mL, were tested in 5 replicates. A 1:100 dilution was performed with one part sample and 99 parts Aptima Specimen Diluent just prior to testing. Testing was performed using one lot of assay reagents on two Panther systems with two Aptima Specimen Diluent lots. The difference between the average reported concentration in native matrix (dilution factor applied to the diluted sample result) and the average concentration in Aptima Specimen Diluent was calculated for each sample set as shown in Table 16 for plasma and Table 17 for serum.

Table 16: Plasma Specimen 1:100 Dilution Matrix Comparison Summary

Plasma Matrix Average Reported Concentration (Log ₁₀ IU/mL) n = 5	Diluent Average Reported Concentration (Log ₁₀ IU/mL) n = 10	Difference Diluent from Plasma Matrix (Log ₁₀ IU/mL)
7.86	7.85	-0.01
7.84	7.83	-0.01
7.78	7.75	-0.03
7.80	7.80	0.00
6.58	6.53	-0.05
6.58	6.52	-0.06
6.58	6.53	-0.05
6.58	6.53	-0.05
9.24 ^a	9.05 ^a	-0.19
9.21 ^a	9.05 ^a	-0.16
9.25 ^a	9.03 ^a	-0.22
9.27 ^a	9.04 ^a	-0.23
9.13ª	8.82ª	-0.31
9.12 ^a	8.81 ^a	-0.31
9.09 ^a	8.84 ^a	-0.25
9.05 ^a	8.84 ^a	-0.21

^aspiked using plasmid DNA

^bn=42

-0.21

-0.20

-0.23

-0.26

-0.29

-0.24

Serum Matrix Diluent Difference **Average Reported Concentration Average Reported Concentration Diluent from Serum Matrix** $(Log_{10} IU/mL) n = 10$ $(Log_{10} IU/mL) n = 5$ (Log₁₀ IU/mL) 7.70 7.85 0.15 7 84 7 85 0.01 7.79 7 82 0.03 7.75 7.79 0.04 6.77 6.77 0.00 6.75 6.80 0.05 6.75 6.71 -0.04 6.70 6.73 0.03 -0 19 9.27^a 9.08^a -0.18 9.24^a 9.06^a

9 08a

9.11^a 8.91^a

8.92^a

8.90^a

8.84^a

Table 17: Serum Specimen 1:100 Dilution Matrix Comparison Summary

9 29^a

9.31^a

9.14^a

9.18^a

9.19^a

9.08^a

Confirmation of the LLoQ in Specimens Diluted in Aptima Specimen Diluent

The LLoQ of the Aptima HBV Quant Assay was confirmed with HBV genotype A clinical specimens diluted into Aptima Specimen Diluent. Specimens were prepared in HBV negative human plasma and serum at 21, 30, and 45 IU/mL. Each panel was diluted 1:3 into Aptima Specimen Diluent just prior to testing to give final concentrations of approximately 7, 10, and 15 IU/mL. A total of 36 replicates of each panel member were tested with one reagent lot across three days. An LLoQ ≤ 10 IU/mL for HBV plasma and serum diluted into Aptima Specimen Diluent was confirmed as shown in Table 18.

Table 18: Confirmation of LLoQ - Samples in Aptima Specimen Diluent

Matrix	% Detected	Aptima HBV Quant (IU/ml)	Aptima HBV Quant (Log ₁₀ IU/ml)	SD (Log ₁₀ IU/ml)	Bias (Log ₁₀ IU/ml)	Calculated TE (Log ₁₀ IU/ml)	Calculated TAE (Log ₁₀ IU/ml)
Plasma	100%	3	0.50	0.19	0.10	0.54	0.48
Serum	100%	2	0.38	0.12	0.46	0.33	0.70

Precision of Diluted Samples

The Aptima HBV Quant precision panel was built by diluting HBV-positive plasma and HBV plasmid DNA into HBV-negative clinical plasma and serum. Positive panels were diluted into Aptima Specimen Diluent. These were tested in five replicates per run by one operator, using three lots of Aptima Specimen Diluent on one Panther system over three test days, two runs a day.

Table 19 shows the precision of assay results (in SD $\log_{10} IU/mL$) for three lots of Aptima Specimen Diluent. Total variability was ≤ 0.15 SD across all panel members and diluent lots.

^aSpiked using plasmid DNA

Aptima® Nonclinical Performance

Table 19: Precision of Panels Diluted in Aptima Specimen Diluent

			Lot 1 Specimen Dili	uent (n=10)	Lot 2 Specimen Dili	ıent (n=10)	Lot 3 Specimen Di	luent (n=10)	Combined Lot	s (n=30)
Matrix	Target Concentration Log ₁₀ IU/mL	Dilution	Average Log ₁₀ IU/mL	SD	Average Log ₁₀ IU/mL	SD	Average Log ₁₀ IU/mL	SD	Average Log ₁₀ IU/mL	SD
	2.20	Neat	3.46	0.07	3.43	0.08	3.46	0.06	3.45	0.07
	3.30	1:3	3.36	0.09	3.35	0.07	3.39	0.09	3.37	0.08
Plasma	4.00	Neat	4.33	0.06	4.27	0.03	4.41	0.05	4.34	0.08
	4.30	1:3	4.34	0.05	4.35	0.05	4.38	0.10	4.35	0.07
	9.18	Neat	9.13	0.05	9.10	0.03	9.26 ^a	0.15	9.16 ^a	0.11
		1:100	9.18	0.03	9.14	0.04	9.33	0.10	9.21	0.10
	0.00	Neat	3.52	0.05	3.48	0.06	3.50	0.07	3.50	0.06
	3.30	1:3	3.45	0.08	3.40	0.06	3.39	0.08	3.41	0.07
Serum	4.00	Neat	4.35	0.05	4.37	0.06	4.43	0.06	4.38	0.06
	4.30	1:3	4.35	0.05	4.37	0.05	4.41	0.04	4.37	0.05
	0.18	Neat	9.08	0.03	9.14	0.05	9.31 ^b	0.12	9.17 ^b	0.12
	9.18 _	1:100	9.18	0.02	9.14	0.03	9.33	0.09	9.22	0.10

SD = standard deviation.

Carryover

To establish that the Panther system minimizes the risk of false positive results arising from carryover contamination, a study was conducted using spiked panels on three Panther systems. Carryover was assessed using high titer HBV DNA spiked plasma samples (8 log₁₀ IU/mL) interspersed between HBV negative samples in a checkerboard pattern. Testing was carried out over fifteen runs. The overall carryover rate was 0.14% (1/705).

^a1 replicate excluded (above calculable range).

^b2 replicates excluded (above calculable range).

Clinical Performance

Clinical Utility

The study was designed to assess the Aptima HBV Quant assay's ability to predict virologic, biochemical, and serologic clinical utility endpoints 48 weeks after initiation of therapy. Specimens were prospectively collected from subjects with chronic HBV infection who were initiating monotherapy with entecavir or tenofovir as part of their standard of care treatment.

Of the 331 subjects enrolled, 86 subjects were not evaluable due to withdrawal, discontinuation, early treatment halt, missing Week 48 results, or missing or low baseline HBV DNA viral loads. The remaining 245 subjects were evaluable for at least one of the clinical utility endpoints and included 126 HBeAg(+) and 119 HBeAg(-) subjects. Table 20 shows the demographic and baseline clinical characteristics of evaluable subjects. The demographic distribution of subjects in this study was consistent with that of patients with chronic HBV in the US.¹⁴ Data from HBeAg(+) and HBeAg(-) subjects were analyzed separately.

Table 20: Demographics and Baseline Clinical Characteristics of Evaluable Subjects

Characteristics		Total
Total, N	N	245
Entecavir	n (%)	94 (38.4)
Tenofovir	n (%)	151 (61.6)
Sex, n (%)	Male	154 (62.9)
	Female	91 (37.1)
Age (years)	Mean ± SD	43.5 ± 13.63
	Median	44.0
	Range	18 – 83
Age category (years), n (%)	18 – 29	40 (16.3)
	30 – 49, n (%)	120 (49.0)
	50 – 70, n (%)	80 (32.7)
	>70, n (%)	5 (2.0)
Ethnicity, n (%)	Hispanic or Latino	7 (2.9)
	Not Hispanic or Latino	236 (96.3)
	Unknown/Refused	2 (0.8)
Race ^a , n (%)	White	89 (36.3)
	Black or African American	16 (6.5)
	Asian	132 (53.9)
	American Indian/Alaska native	0 (0.0)
	Native Hawaiian/Pacific Islander	6 (2.4)
	Other	1 (0.4)
	Unknown/Refused	1 (0.4)
Genotype, n (%)	A	28 (11.4)
	В	64 (26.1)
	С	36 (14.7)
	D	47 (19.2)
	E	2 (0.8)
	F	0 (0.0)
	G	0 (0.0)
	Н	3 (1.2)
	Unknown	65 (26.5)
HBV treatment status, n (%)	Experienced	27 (11.0)
	Naïve	218 (89.0)

Aptima HBV Quant Assay 35 AW-15644, Rev. 001

Clinical Performance Aptima®

Table 20: Demographics and Baseline Clinical Characteristics of Evaluable Subjects (continued)

Characteristics		Total
Previous drug treatment, n (%)	Tenofovir	5 (18.5)
	Entecavir	4 (14.8)
	Adefovir	2 (7.4)
	Lamivudine	1 (3.7)
	Telbivudine	0 (0.0)
	Interferon	10 (37.0)
	Other b	5 (18.5)
Previous treatment outcome, n (%)	Failure	1 (3.7)
	Success	0 (0.0)
	Discontinued for other reasons	26 (96.3)
HBsAg serostatus, n (%)	Positive/reactive	210 (85.7)
	Not done	35 (14.3)
Cirrhotic status, n (%)	Cirrhotic	26 (10.6)
	Non-cirrhotic	201 (82.0)
	Not done	18 (7.3)
HBV viral load (log ₁₀ IU/mL)	Mean ± SD	6.3 ± 1.93
	Median	6.4
	Range	3 – 9
ALT (U/L)	Mean ± SD	102.4 ± 175.97
	Number above ULN ^c	177 (85.9)

^a Subjects may report multiple races.

Prediction of Response to Antiviral Therapy

Aptima HBV Quant assay clinical utility has been assessed for individuals treated with tenofovir and entecavir. No information is available on the assay's clinical utility when other HBV antiviral therapies are used.

Definitions:

Early virologic response outcomes

Week 12 and Week 24 virologic response = HBV DNA <10 IU/mL (<LLoQ) as assessed by the Aptima HBV Quant assay on the Panther system

Week 12 alternative virologic response = HBV DNA ≥2 log₁₀ decrease from baseline

Week 24 alternative virologic response = HBV DNA <2000 IU/mL (for HBeAg+) or <50 IU/mL (for HBeAg-)

Clinical utility endpoints

Week 48 virologic response = HBV DNA <10 IU/mL (<LLoQ) as assessed by an FDA-approved HBV quantitative assay

Alternative Week 48 virologic response = HBV DNA <50 IU/mL as assessed by an FDA-approved HBV quantitative assay

Biochemical response = Normalization of ALT test results at Week 48 (ALT <30 U/L for males and <19 U/L for females)

Serologic response = Loss of HBeAg (HBeAg-negative results) at Week 48

^b Various combinations of the specific drugs listed.

^c The upper limit of normal range (ULN) for alanine aminotransferase (ALT) was 30 U/L for males and 19 U/L for females.

Measures of association and predictive value

Positive Predictive Value (PPV) = True Positive / (True Positive + False Positive) or the probability of response at Week 48 (for the clinical utility endpoint being assessed) in subjects with virologic response at the early time point

Negative Predictive Value (NPV) = True Negative / (False Negative + True Negative) or the probability of non-response at Week 48 (for the clinical utility endpoint being assessed) in subjects with virologic non-response at the early time point

Odds Ratio (OR) = (True Positive × True Negative) / (False Positive × False Negative)

Predicting Week 48 Virologic Response, Defined as HBV DNA <10 IU/mL

In this study, the primary definition of virologic response was HBV DNA <10 IU/mL, and this definition was used for both early virologic response at Weeks 12 and 24, as well as the virologic response at Week 48. The association between early virologic responses at Weeks 12 and 24 and Week 48 clinical utility endpoints (virologic response, biochemical response, serologic response) were assessed.

Predicting Virologic Response at Week 48

Associations between virologic response at Week 48 and virologic response at Week 12 and Week 24 are summarized in Table 21.

Early virologic responses at Weeks 12 and 24, as predictors of Week 48 virologic response, varied by week and treatment.

Table 21: PPV, NPV, and Odds Ratio for Virologic Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <10 IU/mL

			PPV	(%)	NPV	(%)	OR
HBeAg Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^{a,b}
		Entecavir	0.0 (0.0, 93.2)	0/1	82.5 (80.0, 88.4)	33/40	1.49 (<0.01, 30.95)
HBeAg(+) —	12	Tenofovir	100 (27.3, 100)	2/2	74.4 (72.6, 78.7)	61/82	14.30 (1.11, >999.99)
		All	66.7 (15.4, 98.2)	2/3	77.0 (75.7, 79.9)	94/122	6.71 (0.62, 147.55)
	24	Entecavir	50.0 (6.4, 93.2)	2/4	88.2 (83.5, 95.4)	30/34	7.50 (0.74, 79.76)
		Tenofovir	75.0 (52.7, 92.3)	12/16	84.1 (78.5, 90.0)	58/69	15.81 (4.62, 65.54)
		All	70.0 (50.3, 88.1)	14/20	85.4 (81.3, 90.1)	88/103	13.69 (4.74, 44.16
		Entecavir	94.1 (87.1, 99.0)	32/34	22.2 (7.6, 36.0)	4/18	4.57 (0.80, 35.85)
	12	Tenofovir	83.3 (70.5, 93.2)	25/30	46.9 (35.0, 58.9)	15/32	4.41 (1.42, 15.71)
HBeAg(-)		All	89.1 (82.1, 94.7)	57/64	38.0 (29.0, 47.1)	19/50	4.99 (1.96, 14.00)
HBeAg(-) —		Entecavir	93.0 (88.0, 98.1)	40/43	37.5 (6.4, 67.2)	3/8	8.00 (1.21, 55.54)
	24	Tenofovir	82.6 (74.3, 90.4)	38/46	75.0 (54.1, 92.0)	12/16	14.25 (3.92, 62.71)
		All	87.6 (82.7, 92.6)	78/89	62.5 (44.8, 78.0)	15/24	11.82 (4.30, 34.94)

CI=95% profile-likelihood confidence interval

Aptima HBV Quant Assay 37 AW-15644, Rev. 001

^a Shading indicates statistical significance of odds ratios.

^b For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero.

Clinical Performance Aptima®

Predicting Biochemical Response at Week 48

Associations between biochemical response at Week 48 and virologic response at Week 12 and Week 24 are summarized in Table 22.

The value of early virologic responses at Weeks 12 and 24 as a predictor of Week 48 biochemical response varied by week and treatment.

Table 22: PPV, NPV, and Odds Ratio for Biochemical Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <10 IU/mL

			PPV ((%)	NPV ((%)	OR
HBeAg Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^a
	12	Entecavir	NC (NC)	0/0	67.9 (63.4-76.1)	19/28	2.05 (0.01-393.52)
		Tenofovir	100 (27.6, 100)	2/2	58.5 (55.1, 63.9)	31/53	7.00 (0.54, 983.90)
HBeAg(+)		All	100 (27.3, 100)	2/2	61.7 (59.7, 65.5)	50/81	8.02 (0.63, >999.99
	24	Entecavir	66.7 (16.4-98.2)	2/3	68.2 (60.1-80.0)	15/22	4.29 (0.35-101.82)
		Tenofovir	58.3 (33.2, 81.0)	7/12	61.4 (54.3, 69.7)	27/44	2.22 (0.61, 8.61)
		All	60.0 (36.6, 80.9)	9/15	63.6 (58.2, 70.0)	42/66	2.62 (0.84, 8.70)
		Entecavir	43.8 (25.1, 61.2)	7/16	50.0 (25.1, 74.9	6/12	0.78 (0.17, 3.53)
	12	Tenofovir	52.9 (35.1, 72.2)	9/17	76.2 (61.2, 89.8)	16/21	3.60 (0.93, 15.39)
HPoAg()		All	48.5 (36.0, 60.9)	16/33	66.7 (54.5, 78.6)	22/33	1.88 (0.70, 5.20)
HBeAg(-)		Entecavir	45.8 (36.0, 56.0)	11/24	75.0 (25.2, 98.7)	3/4	2.54 (0.28, 55.45)
	24	Tenofovir	42.9 (32.8, 53.1)	12/28	80.0 (53.0, 97.1)	8/10	3.00 (0.61, 22.34)
		All	44.2 (37.8, 50.8)	23/52	78.6 (55.1, 95.0)	11/14	2.91 (0.80, 13.98)

CI=95% profile-likelihood confidence interval, NC=not calculable

Aptima HBV Quant Assay 38 AW-15644, Rev. 001

^a For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero,

Predicting Serologic Response at Week 48

Associations between serologic response at Week 48 and virologic response at Week 12 and Week 24 are summarized in Table 23.

The value of early virologic responses at Weeks 12 and 24 as a predictor of Week 48 serologic response varied by week and treatment.

Table 23: PPV, NPV, and Odds Ratio for Serologic Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <10 IU/mL

			PPV	(%)	NPV	(%)	OR
HBeAq Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^a
	12	Entecavir	100 (6.5, 100)	1/1	86.8 (84.2, 93.9)	33/38	18.27 (0.86, >999.99)
		Tenofovir	0.0 (0.0, 72.0)	0/2	82.9 (81.7, 86.0)	68/82	0.95 (<0.01, 12.45)
HBeAg(+) _		All	33.3 (1.8, 84.3)	1/3	84.2 (83.1, 86.8)	101/120	2.66 (0.12, 29.11)
пведу(т) –		Entecavir	50.0 (6.4, 93.2)	2/4	88.2 (83.5, 95.4)	30/34	7.50 (0.74, 79.76)
	24	Tenofovir	18.8 (3.1, 39.1)	3/16	84.1 (80.6, 89.1)	58/69	1.22 (0.25, 4.59)
		All	25.0 (8.5, 43.6)	5/20	85.4 (82.5, 89.4)	88/103	1.96 (0.57, 5.94)

CI=95% profile-likelihood confidence interval

Aptima HBV Quant Assay 39 AW-15644, Rev. 001

^a For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero.

Clinical Performance Aptima®

Predicting Week 48 Virologic Response, Defined as HBV DNA <50 IU/mL (Alternative Definition)

Alternative definitions of early (Weeks 12 and 24) and Week 48 virologic responses also were assessed (see alternative response definitions above).

Associations between clinical utility endpoints and virologic response at Week 12 and Week 24, using these alternate definitions of virologic response, are summarized in Table 24 (virologic response), Table 25 (biochemical response), and Table 26 (serologic response).

Table 24: PPV, NPV, and Odds Ratio for Virologic Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <50 IU/mL

			PPV	(%)	NPV	(%)	OR
HBeAg Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^{a,b}
		Entecavir	34.2 (26.9, 39.2)	13/38	66.7 (16.1, 98.2)	2/3	1.04 (0.09, 23.60)
	12	Tenofovir	55.1 (51.9, 59.4)	43/78	83.3 (43.5, 99.2)	5/6	6.14 (0.93, 120.54)
LIDo Ag(1)		All	48.3 (45.6, 51.3)	56/116	77.8 (45.3, 97.1)	7/9	3.27 (0.75, 22.54)
HBeAg(+) —	24	Entecavir	65.0 (51.3, 81.2)	13/20	100 (86.3, 100)	18/18	66.59 (7.20, >999.99)
		Tenofovir	72.4 (64.7, 80.6)	42/58	88.9 (74.6, 97.2)	24/27	21.00 (6.28, 97.36)
		All	70.5 (63.5, 77.9)	55/78	93.3 (84.0, 98.3)	42/45	33.47 (10.81, 148.13)
		Entecavir	100 (NC)	52/52	NC (NC)	0/0	NC
	12	Tenofovir	93.0 (89.1, 97.5)	53/57	60.0 (21.3, 93.3)	3/5	19.87 (2.62, 191.49)
HReAg(-)		All	96.3 (94.4, 98.7)	105/109	60.0 (21.1, 93.3)	3/5	39.37 (5.24, 376.77)
HBeAg(-) —		Entecavir	100 (NC)	47/47	0.0 (NC)	0/4	NC
	24	Tenofovir	93.9 (88.6, 98.3)	46/49	30.8 (6.7, 52.4)	4/13	6.81 (1.30, 39.97)
		All	96.9 (93.8, 99.2)	93/96	23.5 (7.0, 39.8)	4/17)	9.54 (1.91, 53.22)

CI=95% profile-likelihood confidence interval, NC=not calculable

^a Shading indicates statistical significance of odds ratios.

^b For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero, unless there were either no Week 48 responses or no Week 48 non-responses, which resulted in reporting the odds ratio as NC.

Aptima® Clinical Performance

Table 25: PPV, NPV, and Odds Ratio for Biochemical Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <50 IU/mL

			PP	V	NP	/	OR
HBeAg Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^{a,b}
		Entecavir	33.3 (25.1-39.0)	9/27	100 (6.7-100)	1/1	1.54 (0.07-233.77)
	12	Tenofovir	44.2 (39.6, 48.6)	23/52	66.7 (15.7, 98.2)	2/3	1.59 (0.14, 35.36)
HBeAg(+) –		All	40.5 (37.1, 43.7)	32/79	75.0 (23.9, 98.7)	3/4	2.04 (0.25, 42.30)
	24	Entecavir	50.0 (31.2-69.5)	7/14	81.8 (58.9-97.1)	9/11	4.50 (0.79-37.15)
		Tenofovir	52.5 (44.0, 61.8)	21/40	81.3 (60.1-96.7)	13/16	4.79 (1.30, 23.28)
		All	51.9 (44.1, 60.2)	28/54	81.5 (66.4, 93.1)	22/27	4.74 (1.66, 15.82)
		Entecavir	46.4 (39.5, 52.6)	13/28	NC (NC)	0/0	0.87 (<0.01, 166.17)
	12	Tenofovir	40.0 (33.6, 46.3)	14/35	100 (41.4, 100)	3/3	4.72 (0.41, 653.11)
UPoAg()		All	42.9 (39.6, 46.7)	27/63	100 (41.2, 100)	3/3	5.27 (0.48, 720.38)
HBeAg(-) —		Entecavir	44.0 (34.4, 52.7)	11/25	66.7 (16.3-98.2)	2/3	1.57 (0.13, 36.42)
	24	Tenofovir	41.4 (31.3, 51.1)	12/29	77.8 (48.6, 97.1)	7/9	2.47 (0.49, 18.57)
		All	42.6 (36.4, 48.9)	23/54	75.0 (48.8, 93.7)	9/12	2.23 (0.59, 10.87)

CI=95% profile-likelihood confidence interval

Table 26: PPV, NPV, and Odds Ratio for Serologic Response Predicted by Early Virologic Response During Treatment: Week 48 Virologic Response Defined as <50 IU/mL

			PP	V	NP'	V	OR
HBeAg Status	Week of Early Virologic Response	Treatment	Estimate (95% CI)	n/N	Estimate (95% CI)	n/N	Estimate (95% CI) ^{a,b}
	12	Entecavir	16.7 (10.4, 19.6)	6/36	100 (44.1, 100)	3/3	1.49 (0.12, 209.89)
		Tenofovir	16.7 (12.9, 18.6)	13/78	83.3 (44.2, 99.2)	5/6	1.00 (0.14, 19.98)
HBeAg(+)		All	16.7 (13.9, 18.1)	19/114	88.9 (58.8, 99.5)	8/9	1.60 (0.27, 30.56)
ПВеду(+)		Entecavir	30.0 (16.5, 41.6)	6/20	100 (86.3, 100)	18/18	16.59 (1.72, >999.99)
	24	Tenofovir	22.4 (16.9, 26.8)	13/58	96.3 (84.7, 99.9)	26/27	7.51 (1.37, 140.31)
		All	24.4 (20.0, 28.4)	19/78	97.8 (90.4, 99.9)	44/45	14.17 (2.77, 259.25)

CI=95% profile-likelihood confidence interval

Aptima HBV Quant Assay 41 AW-15644, Rev. 001

^a Shading indicates statistical significance of odds ratios.

^b For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero.

^a Shading indicates statistical significance of odds ratios.

^b For calculation of odds ratios and their confidence intervals, 0.5 was added to all cells whenever at least one cell was zero.

Clinical Performance Aptima®

Conclusion

Overall, the results demonstrate that the Aptima HBV Quant assay can quantitate HBV DNA levels at baseline and during treatment to aid in assessing viral response to treatment. This study demonstrated that early virologic responses at Weeks 12 and 24, as predictors of Week 48 virologic response, varied by week and treatment.

The Aptima HBV Quant assay can be used as an aid in the management of chronic HBV-infected patients undergoing HBV antiviral drug therapy.

Reproducibility

Reproducibility was evaluated on the Panther system at three external U.S. sites. Two operators performed testing at each site. Each operator performed two runs per day over three days, using three reagent lots over the course of testing. Each run had three replicates of each panel member. Overall, 108 replicates of each panel member were tested.

Reproducibility was tested using panel members prepared using HBV-negative plasma. The positive panel members were positive for HBV genotypes A or C. HBV DNA concentrations spanned the linear range of the assay.

Table 27 shows the reproducibility and precision of assay results for each positive panel member between sites, between operators/days, between lots, between runs, within runs, and overall.

The coefficient of variation was calculated using the following equation where σ^2 is the sample variance of the data after \log_{10} transformation.

$$%CV = 100 \times \sqrt{(10^{\sigma^2 \ln(10)} - 1)}$$

For all HBV-positive panel members, agreement values were 100%.

Aptima® Clinical Performance

Table 27: Reproducibility of Aptima HBV Quant Assay HBV DNA Levels on the Panther System in Positive Panel Members

		Obse	rved Mean			oution to Total Varia SD (%CV)	nce					
GT	N	IU/mL	Log ₁₀ IU/mL	Between Sites	Between Operators/Days ^a	Between Lots	Between Runs	Within Runs	Total Variance SD (%CV)			
	108	17.6	1.2	0.059	<0.001	0.138	0.090	0.178	0.250			
	100	17.0	1.2	(13.578)	(<0.001)	(32.693)	(20.869)	(42.883)	(62.666)			
	108	120.4	2.1	0.009	0	0.074	0.051	0.106	0.139			
	100	129.4	2.1	(2.162)	(0)	(17.109)	(11.869)	(24.736)	(32.886)			
	107	1056.0	3.0	0.035	0.032	0.014	0.032	0.085	0.103			
	107	0.0001	3.0	(7.994)	(7.432)	(3.246)	(7.356)	(19.666)	(24.060)			
	100	7663.0	2.0	0	0.027	0.040	0.044	0.066	0.092			
Α	108	7003.0	3.9	(0)	(6.262)	(9.235)	(10.088)	(15.194)	(21.540)			
	108	100170 1	E 2	0.027	0.042	0.042	0.030	0.072	0.102			
	106	188172.1	5.3	(6.281)	(9.707)	(9.787)	(6.829)	(16.689)	(23.772)			
	400	0200004.4	7.0	0.038	0	0.031	0.064	0.068	0.106			
	108	9389094.1	7.0	(8.846)	(0)	(7.237)	(14.791)	(15.756)	(24.692)			
	107	06664677.0	7.0	0.038	0.029	0.020	0.037	0.049	0.081			
	107	86664677.2	7.9	(8.692)	(6.753)	(4.584)	(8.635)	(11.375)	(18.725)			
	407	750700400.0	0.0	0.024	0.052	0.015	0.045	0.053	0.091			
	107	753726183.2	8.9	(5.476)	(11.997)	(3.499)	(10.304)	(12.187)	(21.163)			
	107	17.0	17.0	17.0	17.0	1.0	0.041	0.041	0.074	0.092	0.189	0.230
	107	17.0	17.0 1.2	(9.521)	(9.392)	(17.147)	(21.438)	(45.704)	(57.010)			
	400	450.0	0.0	0.035	0	0.055	0.064	0.131	0.160			
	108	152.9	2.2	(8.127)	(0)	(12.706)	(14.925)	(30.748)	(38.013)			
	108	1363.8	3.1	0.042	0.023	0	0.061	0.055	0.094			
	100	1303.0	3.1	(9.583)	(5.316)	(0)	(14.033)	(12.623)	(22.002)			
	108	9871.9	4.0	0.011	0.014	0.040	0.038	0.059	0.083			
С	100	9071.9	4.0	(2.472)	(3.270)	(9.337)	(8.801)	(13.651)	(19.291)			
	400	247400 5	F 2	0.031	0.047	0.016	0.026	0.063	0.090			
	108	217400.5	5.3	(7.255)	(10.843)	(3.791)	(6.023)	(14.685)	(21.044)			
	100	12007170 F	7.4	0.046	0	0.020	0.064	0.073	0.109			
	108	12087179.5	7.1	(10.543)	(0)	(4.652)	(14.762)	(16.922)	(25.501)			
•	100	E7742712 0	7.8	0.044	0.028	0.013	0.043	0.052	0.087			
	۱۵۵	108 57743712.8	1.8	(10.232)	(6.472)	(2.944)	(10.026)	(12.010)	(20.146)			
	400	F704047F4 0	0.7	0.042	0.048	0.028	0.034	0.048	0.091			
	108	572184754.9	8.7	(9.711)	(11.160)	(6.374)	(7.740)	(11.081)	(21.208)			

[%]CV=log-normal coefficient of variation, GT=genotype, SD=standard deviation

Note: Variability from some factors may be numerically negative. This can occur if the variability due to those factors is very small. In these cases, SD and CV are shown as 0.

Aptima HBV Quant Assay 43 AW-15644, Rev. 001

^a Between Operators may be confounded with Between Days, therefore, Between Operators and Between Days estimates are combined in Between Operators/Days.

Bibliography Aptima®

Bibliography

1. **Terrault NA, Bzowej NH, Chang KM, Hwang JP, Jonas MM, Murad MH**. AASLD Guidelines for Treatment of Chronic Hepatitis B. Hepatology. 2016; 63(1):261-283. www.aasld.org/sites/default/files/guideline_documents/hep28156.pdf. Accessed July 2017.

- 2. **Aspinall EJ, Hawkins G, Fraser A, Hutchinson SJ, Goldberg D.** Hepatitis B prevention, diagnosis, treatment and care: a review. Occupational Medicine 2011; 61(8):531-540.
- 3. **Centers for Disease Control and Prevention/National Institutes of Health.** Early Identification and Linkage to Care of Persons with Chronic Hepatitis B Virus Infection Three U.S. Sites, 2012-2014. May 9, 2014, 63(18);399-401.
- 4. Liaw YF, Chu CM. Hepatitis B virus infection. Lancet. 2009;373(9663):582-592.
- 5. **European Association for the Study of the Liver.** EASL Clinical Practice Guidelines: Management of chronic hepatitis B virus infection. Journal of Hepatology 2012; 57:167-185.
- 6. International Agency for Research on Cancer. Hepatitis B Virus. IARC Monographs 2012; 100B: 93-133.
- 7. Price, J. An Update on Hepatitis B, D, and E Viruses. Topics in Antiviral Medicine 2014; 21(5): 157-163.
- 8. Clinical and Laboratory Standards Institute. 2005. Collection, Transport, Preparation, and Storage of Specimens for Molecular Methods; Approved Guideline. CLSI Document MM13-A. Wayne, PA.
- 9. 29 CFR Part 1910.1030. Occupational Exposure to Bloodborne Pathogens; current version.
- Centers for Disease Control and Prevention/National Institutes of Health. Biosafety in Microbiological and Biomedical Laboratories (BMBL); current version.
- 11. Clinical and Laboratory Standards Institute. 2002. Clinical Laboratory Waste Management. CLSI Document GP5-A2. Villanova, PA.
- 12. Clinical and Laboratory Standards Institute (CLSI). 2012. Evaluation of Detection Capability for Clinical Laboratory Measurement Procedures; Approved Guideline—Second Edition. CLSI Document EP17-A2. Clinical and Laboratory Standards Institute, Wayne, PA.
- 13. Clinical and Laboratory Standards Institute (CLSI). 2003. Evaluation of the Linearity of Quantitative Measurement Procedures: A Statistical Approach; Approved Guideline. CLSI document EP06-A. Clinical and Laboratory Standards Institute, Wayne, PA
- Ghany MG, Perrillo R, Li R, et al. Characteristics of adults in the Hepatitis B Research Network in North America reflect their country of origin and hepatitis B virus genotype. Clin Gastroenterol Hepatol. 2015;13(1):183-192. doi:http://dx.doi.org/10.1016/j.cgh.2014.06.028.





Hologic, Inc. 10210 Genetic Center Drive San Diego, CA 92121 USA

Customer Support: +1 800 442 9892

customersupport@hologic.com

Technical Support: +1 888 484 4747

molecularsupport@hologic.com

For more contact information, visit www.hologic.com.

Hologic, Aptima, Panther, and associated logos are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries. All other trademarks that may appear in this package insert are the property of their respective owners.

This product may be covered by one or more U.S. patents identified at www.hologic.com/patents.

© 2017-2018 Hologic, Inc. All rights reserved. AW-15644, Rev. 001 2018-01