



Definity™

Cervical Dilator

with SureAccess™ Balloon Technology




REDEFINE DILATION

Redefine Dilation

Definity cervical dilator with **SureAccess** balloon technology

Effortlessly access a wide variety of cavities with innovative balloon technology that offers a gentle insertion and dilation you can be confident in.

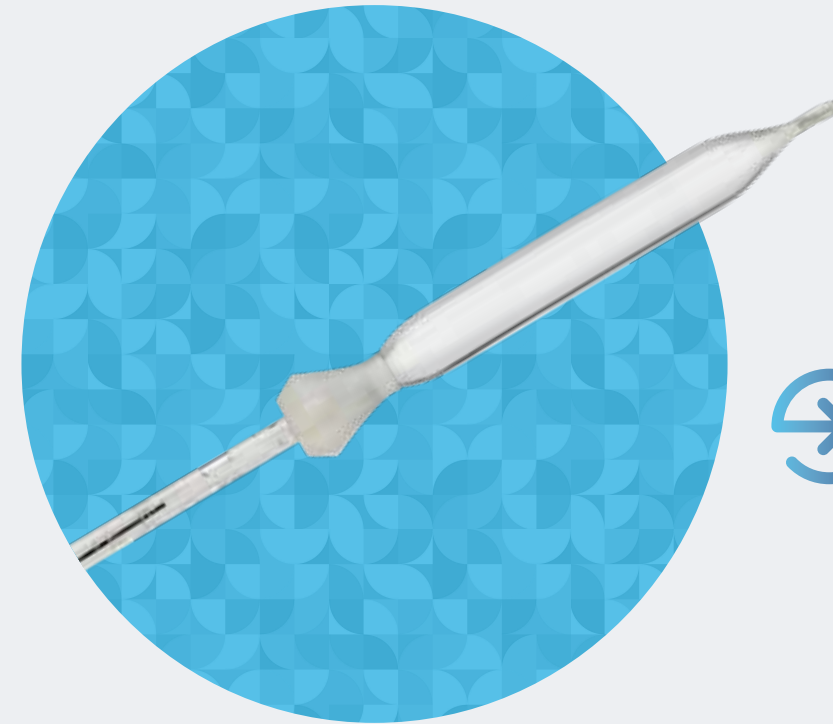

VISUAL INDICATORS

Color coding allows the user to view progress & charging of the device

DILATION LEVER

Easily control when dilation begins

ACORN SEAL

Enables secure positioning

SUREACCESS BALLOON TECHNOLOGY

Designed to dilate even stenotic or challenging anatomies¹



INFLATION DEVICE

Reduces hand fatigue by maintaining pressure for the physician

ERGONOMIC DESIGN

Comfortable to hold for left or right handed users

SEAMLESS Expandable balloon designed to eliminate multiple passes

GENTLE No tenaculum needed for reduced cervical trauma

CONSISTENT Enable dilation through challenging anatomies, including stenosis¹

CONFIDENT Reduce risk of perforation and creation of false passage²

Experience the next major advancement in dilation in over 140 years!

3 sizes to meet your hysteroscopic needs, enabling dilation through challenging anatomies, including stenosis, tortuous canals or cervixes associated with severely anteverted or retroverted uteri.¹

Below are a few possible opportunities in which the Definity cervical dilator can be used. Please note that this list is not exhaustive.

5mm	7mm	9mm
IUD Placement	NovaSure® Advanced device	NovaSure® Classic device
Endometrial Sampling	Omni™ hysteroscope	Omni™ hysteroscope
	MyoSure® Classic hysteroscope	MyoSure® XL hysteroscope

ORDERING INFORMATION

Description	Part #
5mm Definity cervical dilator, 5 pack	DCD 505
7mm Definity cervical dilator, 5 pack	DCD 705
9mm Definity cervical dilator, 5 pack	DCD 905

Experience innovative dilation for yourself, contact your rep or visit: gynsurgicalsolutions.com/definity

REFERENCES

1. Hologic Data on File. VER-09559. Bench test utilizing 2.4mm model.
2. As compared to fixed dilators.

IMPORTANT SAFETY INFORMATION

The Definity™ cervical dilator catheter system is intended to be used whenever cervical softening and dilation is desired. Some examples are: treatment of cervical stenosis, IUD placement and removal, placement of instruments for intrauterine radiotherapy, endometrial biopsy, global endometrial ablation, uterine tissue removal, uterine curettage, diagnostic hysteroscopy, operative hysteroscopy. This device is not intended for use in the induction of labor. Use of the Definity cervical dilator catheter system is contraindicated in patients with: an active genital tract infection such as genital herpes, pelvic structure abnormality that prevents passage of the device, or invasive cervical cancer. This device is also contraindicated for the induction of labor.